

# Phantom SDK

## Reference Manual

---

Version 13.4.788.0 (September 2019)

Copyright (C) 1992-2019 Vision Research Inc. All Rights Reserved.

*Notice*

The information contained in this document file includes data that is proprietary of Vision Research Inc. and shall not be duplicated, used, or disclosed – in whole or in part – without the prior written permission of Vision Research Inc. The data subject to this restriction is contained in all pages of this file.

THIS PUBLICATION AND THE INFORMATION HEREIN ARE FURNISHED AS IS, ARE FURNISHED FOR INFORMATIONAL USE ONLY, ARE SUBJECT TO CHANGE WITHOUT NOTICE, AND SHOULD NOT BE CONSTRUED AS A COMMITMENT BY VISION RESEARCH INC. VISION RESEARCH INC ASSUMES NO RESPONSIBILITY OR LIABILITY FOR ANY ERRORS OR INACCURACIES THAT MAY APPEAR IN THE INFORMATIONAL CONTENT CONTAINED IN THIS GUIDE, MAKES NO WARRANTY OF ANY KIND (EXPRESS, IMPLIED, OR STATUTORY) WITH RESPECT TO THIS PUBLICATION, AND EXPRESSLY DISCLAIMS ANY AND ALL WARRANTIES OF MERCHANTABILITY, FITNESS FOR PARTICULAR PURPOSES, AND NONINFRINGEMENT OF THIRD-PARTY RIGHTS.

*Licenses and Trademarks*

Windows is a trademark of Microsoft Corporation.

*Software release*

Some new fields of the SETUP structure may be added in new software releases. This document is based on software release 788 (PhCon.dll, PhFile.dll, PhInt.dll version 13.4. 788.0, PCC version 3.4.788.0).



**AMETEK**<sup>®</sup>  
MATERIALS ANALYSIS DIVISION

100 Dey Road  
Wayne, NJ 07470 USA  
+1.973.696.4500  
[phantom@visionresearch.com](mailto:phantom@visionresearch.com)

[www.phantomhighspeed.com](http://www.phantomhighspeed.com)

## Contents

<b>1. Overview .....</b>	<b>8</b>
1.1. About This Document .....	8
1.2. Phantom SDK Structure .....	8
1.3. Backward Compatibility .....	9
1.4. System Requirements .....	9
1.5. Supported Phantom Cameras .....	10
<b>2. Introduction .....</b>	<b>11</b>
2.1. Entities .....	11
2.1.1. The Camera .....	11
2.1.2. The Cine .....	13
2.1.3. Image Numbering .....	13
2.2. Relations .....	14
2.2.1. The Camera Pool .....	14
2.2.2. Cameras and Cines .....	14
2.3. Actions .....	15
2.3.1. Camera and Cine Control .....	15
2.3.2. Cine Use Case .....	16
2.3.3. Putting It All Together - Creating a New Application .....	16
<b>3. Data Structure Reference .....</b>	<b>18</b>
3.1. Time Related Data Structures .....	18
3.1.1. TIME64 .....	18
3.1.2. TC .....	19
3.1.3. TCU .....	20
3.2. Acquisition Parameters .....	21
3.2.1. ACQUIPARAMS .....	21
3.3. Cine Related Data Structures .....	25
3.3.1. IMRANGE .....	25
3.3.2. CINESTATUS .....	25
3.4. Camera Related Data Structures .....	26
3.4.1. CAMERAOPTIONS .....	26
3.4.2. CAMERAID .....	30
3.4.3. PULSEPARAM .....	31
3.5. Save Related Data Structures .....	32
3.5.1. SAVECINEINFO .....	32
3.6. Image Headers .....	32
3.6.1. BITMAPINFOHEADER .....	32
3.6.2. IH .....	34
3.7. Image Processing Related Data Structures .....	35
3.7.1. WBGAIN .....	35
3.7.2. IMFILTER .....	35
3.7.3. REDUCE16TO8PARAMS .....	36
3.7.4. TONEDESC .....	36
3.7.5. CMDESC .....	36
<b>4. User Defined Functions .....</b>	<b>37</b>

4.1. PROGRESSCALLBACK .....	37
4.2. PROGRESSCB .....	37
4.3. DEMOSAICINGFUNCPTR .....	38
<b>5. Function Reference .....</b>	<b>39</b>
5.1. Initialize/Finalize.....	39
5.1.1. PhRegisterClientEx .....	39
5.1.2. PhUnregisterClient.....	39
5.2. Create & Destroy Cine Objects.....	40
5.2.1. PhGetCineLive.....	40
5.2.2. PhNewCineFromCamera .....	40
5.2.3. PhNewCineFromFile .....	41
5.2.4. PhDuplicateCine .....	41
5.2.5. PhGetOpenCineName.....	42
5.2.6. PhDestroyCine .....	42
5.3. Read Cine Images.....	43
5.3.1. PhGetCineImage .....	43
5.3.2. PhSetUseCase.....	44
5.3.3. PhGetUseCase .....	44
5.4. Save Cine Images to File .....	45
5.4.1. PhGetSaveCineName .....	45
5.4.2. PhGetSaveCinePath .....	45
5.4.3. PhWriteCineFile .....	46
5.4.4. PhWriteCineFileAsync.....	46
5.4.5. PhStopWriteCineFileAsync .....	47
5.4.6. PhGetWriteCineFileProgress .....	47
5.4.7. PhGetWriteCineFileProgressAll.....	48
5.5. Cine Parameters .....	48
5.5.1. PhSetCineInfo.....	48
5.5.2. PhGetCineInfo.....	49
5.5.3. PhCineGet .....	49
5.5.4. PhCineSet .....	50
5.5.5. PhGetCineAuxData .....	50
5.5.6. PhMeasureCineWB.....	51
5.5.7. PhPrintTime .....	52
5.5.8. PhGetCineTriggerTime .....	53
5.6. Acquisition Parameters .....	53
5.6.1. PhGetCineParams .....	53
5.6.2. PhSetCineParams .....	54
5.6.3. PhSetSingleCineParams.....	54
5.6.4. PhGetBitDepths.....	55
5.6.5. PhGetResolutions .....	55
5.6.6. PhGetFrameRateRange .....	57
5.6.7. PhGetExactFrameRate.....	57
5.6.8. PhGetExposureRange.....	58
5.6.9. PhParamsChanged.....	58
5.7. Camera Pool.....	59
5.7.1. PhGetCameraCount .....	59
5.7.2. PhConfigPoolUpdate.....	59
5.7.3. PhGetAllIpCameras.....	60
5.7.4. PhOffline .....	60
5.7.5. PhMakeAllIpVisible.....	60

5.7.6.	PhGetVisibleIp .....	61
5.7.7.	PhAddVisibleIp .....	61
5.7.8.	PhRemoveVisibleIp .....	62
5.7.9.	PhDisableVisible .....	62
5.7.10.	PhIsVisibleDisabled .....	62
5.7.11.	PhFillIDInfo .....	63
5.7.12.	PhGetIgnoredIp .....	63
5.7.13.	PhAddIgnoredIp .....	63
5.7.14.	PhRemoveIgnoredIp .....	64
5.7.15.	PhAddSimulatedCamera .....	64
5.8.	Camera Parameters .....	65
5.8.1.	PhGet .....	65
5.8.2.	PhSet .....	65
5.8.3.	PhGet1 .....	65
5.8.4.	PhSet1 .....	66
5.8.5.	PhGet2 .....	67
5.8.6.	PhGetCameraOptions .....	67
5.8.7.	PhSetCameraOptions .....	68
5.8.8.	PhGetVersion .....	68
5.8.9.	PhGetCameraID .....	69
5.8.10.	PhSetCameraName .....	69
5.8.11.	PhFindCameraNumber .....	70
5.8.12.	PhGetCameraModel .....	70
5.8.13.	PhGetCameraTime .....	70
5.8.14.	PhSetCameraTime .....	71
5.8.15.	PhMaxCineCnt .....	71
5.8.16.	PhFirstFlashCine .....	72
5.8.17.	PhGetPartitions .....	72
5.8.18.	PhSetPartitions .....	72
5.8.19.	PhGetCineStatus .....	73
5.8.20.	PhGetCineCount .....	73
5.8.21.	PhMemorySize .....	74
5.9.	Camera Control .....	74
5.9.1.	PhRecordCine .....	74
5.9.2.	PhRecordSpecificCine .....	74
5.9.3.	PhSendSoftwareTrigger .....	75
5.9.4.	PhDeleteCine .....	75
5.9.5.	PhVideoPlay .....	76
5.9.6.	PhGetVideoFrNr .....	76
5.9.7.	PhGetCameraErrMsg .....	77
5.10.	Camera Calibration .....	77
5.10.1.	PhBlackReference .....	77
5.10.2.	PhBlackReferenceCI .....	78
5.10.3.	PhWriteStgFlash .....	79
5.11.	NVM .....	79
5.11.1.	PhNVMGetStatus .....	79
5.11.2.	PhNVMerase .....	80
5.11.3.	PhNVMContRec .....	80
5.11.4.	PhNVMGetSaveRange .....	81
5.11.5.	PhNVMSetSaveRange .....	81
5.11.6.	PhNVMSave .....	82
5.11.7.	PhNVMSaveClip .....	82
5.11.8.	PhNVMRestore .....	83

5.12.	DLLs Options .....	83
5.12.1.	PhSetDllsOption.....	83
5.12.2.	PhGetErrorMessage.....	84
5.12.3.	PhSetDllsLogOption.....	84
5.12.4.	PhGetDllsLogOption .....	85
5.13.	Other - Image Handling & Processing.....	85
5.13.1.	PhIHtoBITMAPINFOHEADER .....	85
5.13.2.	PhProcessImage .....	86
5.13.3.	PhInterpolateColor .....	87
5.13.4.	PhTriplicateGray .....	88
5.13.5.	PhWBAdjustUnint.....	88
5.13.6.	PhComputeWB .....	89
5.13.7.	PhImHisto.....	90
5.14.	LabView Specific Functions .....	91
5.14.1.	PhLVRegisterClientEx .....	91
5.14.2.	PhLVUnregisterClient.....	91
5.14.3.	PhLVProgress .....	92
5.15.	Nucleus Specific Functions.....	94
5.15.1.	PhUpgradeFirmware.....	94
<b>6.</b>	<b>File Naming Convention .....</b>	<b>95</b>
6.1.	Converting Cines into Multiple Images.....	95
6.1.1.	Converting Single Cine into Multiple Images .....	95
6.1.2.	Converting Multiple Cines into Multiple Images.....	96
6.2.	Examples.....	96
<b>7.</b>	<b>Supported File Formats.....</b>	<b>99</b>
<b>8.</b>	<b>Error and Warning Codes .....</b>	<b>100</b>
8.1.	Error Codes.....	100
8.2.	Warning Codes .....	104
<b>Appendix 1.</b>	<b>PhGetCineInfo/PhSetCineInfo Selectors.....</b>	<b>107</b>
<b>Appendix 2.</b>	<b>Image Filtering.....</b>	<b>115</b>
<b>Appendix 3.</b>	<b>PhGet/PhSet Selectors .....</b>	<b>117</b>
<b>Appendix 4.</b>	<b>PhCineGet/PhCineSet Selectors.....</b>	<b>123</b>
<b>Appendix 5.</b>	<b>PhGet1/PhSet1 Selectors.....</b>	<b>125</b>
<b>Appendix 6.</b>	<b>PhGet2 Selectors.....</b>	<b>126</b>
<b>Appendix 7.</b>	<b>PhUpgradeFirmware Selectors .....</b>	<b>127</b>



## 1. Overview

### 1.1. About This Document

This document provides information for software developers who are writing applications using Phantom SDK for controlling Vision Research high-speed cameras and who need to:

- configure cameras and set acquisition parameters
- send record/trigger commands to a camera
- obtain live or recorded images from a camera
- save images in different video file formats

Chapter 2 introduces the Phantom camera general model thoroughly describing each component element, interactions between these elements and possible control operations. It is not intended to be a presentation of Vision Research cameras and their features, but a conceptual description of the software perspective over these cameras' universe. You can skip this chapter if you already are familiar with notions such as *camera*, *cine*, *cinehandle* or *Frame Buffer Memory*.

Chapter 3 describes Phantom SDK data structures, chapter 4 lists user defined function types, while chapter 5 details Phantom SDK functions. Please note that some of the data structures and functions available in Phantom SDK headers may not be documented here. They are deprecated and while best effort is put to maintain correct functionality, their use is not recommended.

### 1.2. Phantom SDK Structure

The Phantom SDK includes:

- this reference manual
- 32-bit and 64-bit binary .dll files
- public header .h files for C/C++
- public header .cs files for C#
- 32-bit and 64-bit static library .lib files
- Demo C/C++
  - HelloPhantomC - a C/C++ hello world demo application
  - PhSimpleDemoCPP - a C/C++ simple GUI demo application
  - PhDemoCPP - a C/C++ GUI demo application
- Demo C#
  - PhSimpleDemoCS - a C# simple GUI demo application
  - PhDemoCS - a C# GUI demo application
- Demo Matlab
  - HelloPhantom - a Matlab hello world demo application
  - SimplePhMatlabDemo - a Matlab simple GUI demo application
  - PhDemoMathlab - a Matlab GUI demo application
  - SimpleCineFileReader - a Matlab read cine file demo application
- Demo LabView
  - SimpleVideoAnalyzer - a LabView demo application
  - Automated Save - a LabView demo to save Cine files



Current SDK version is available for Microsoft operating systems only, both 32 and 64 bit versions. All code has been written and tested under Microsoft Visual Studio 2013. Phantom SDK functions can also be used from Matlab or LabView applications as depicted by the two included demos.

There are six main dynamic link libraries: *PhCon*, *PhFile*, *PhInt*, *PhSig*, *PhSigV* and *PhRange*. Each of them deals with different camera control components. However, they are interdependent and must therefore be used together.

*PhCon* is the control library for Phantom high-speed cameras. It allows setting camera parameters, recording a movie in the camera *Frame Buffer Memory* (FBM) or transferring images to the client's computer through camera's Ethernet connection or high speed 1394 bus.

*PhFile* is the library responsible for files and movies management. It offers support for reading, saving, converting recordings no matter if stored on camera or inside a file.

*PhInt* is the color interpolation and image processing library for Phantom cameras. It implements various algorithms for demosaicing raw sensor information in order to output RGB color images.

*PhSig* and *PhSigV* deal with signal acquisition using a separate and optional device: The Signal Acquisition Module (SAM) produced by Vision Research. This module allows acquiring analog and binary signals, in parallel and synchronized with the image acquisition. It is controlled from Vision Research applications only (PCC or Phantom) therefore these dlls are not described in this document.

*PhRange* deal with signal acquisition on the range data input. This module allows acquiring the range data signals, in parallel and synchronized with the image acquisition.

### 1.3. Backward Compatibility

The current version of the libraries is a superset of previous versions, made public starting in May 2001. It maintains compatibility with old third party applications written for a previous SDK version. Such an application can be updated simply by replacing old binary Phantom files with the latest version. The update enables the application to get access to new cameras previous unavailable. However, in order to fully take advantage of new cameras and new features, a more complex update is required involving, among other steps, a full application recompiling.

### 1.4. System Requirements

It is recommended that the computer you use have a Pentium class microprocessor running 1.7 GHz or higher. The Phantom SDK can be used on Microsoft operating systems only:

- Windows 10 (32 & 64 bit)
- Windows 8/8.1 (32 & 64 bit)
- Windows 7 (32 & 64 bit)
- Windows Vista (32 & 64 bit)

- Windows 2000

A more powerful computer and fast disk drive will deliver faster display and playback speeds as well as shorter save and download times for your application. Also, image processing functions are configured to take advantage of multicore processors.

All recent Vision Research cameras offer a 1Gbit/s Ethernet connection for both camera control and image data transfer; therefore, a corresponding Ethernet card would be required for your computer. For Cine Stations (readers for our CineMags) and some cameras an optional 10Gbit/s Ethernet connection is available enabling even faster transfer rates.

### **1.5. Supported Phantom Cameras**

All Vision Research high speed video cameras produced until September, 2019 except multi-head cameras and I3 external memories are supported by this SDK version. However, please note that the software has not been tested on older cameras, such as firewire cameras.

## 2. Introduction

The Phantom SDK provides a way to programmatically access and control Phantom cameras and recorded video movies. This programming model consists of *entities*, *relations* between entities and *actions*.

### 2.1. Entities

There are two main types of entities the Phantom SDK is working with: *camera* entities and *cine* entities.

#### 2.1.1. The Camera

A *camera* is a generalization of an actual Phantom camera. A *camera* allows a unitary, simplified approach no matter if working with the latest generation Vision Research camera or with other older models.

A camera is responsible for recording and storing images. It performs these actions following explicit requests and conforming to pre-configured parameters. There are four main types of parameters:

- camera options – used for general camera configuration (for example: what video output standard to use, should the camera store acquired images on the CineMag if one is present, what should the camera do after a recording finishes)
- acquisition parameters – controlling a movie acquisition (for example: frame rate, exposure time, image resolution)
- image processing parameters – various image adjustment parameters (for example: brightness, saturation, hue)
- other parameters – parameters not necessarily in one of the above categories (for example: lens information)

The first two types of parameters are represented by the `CAMERAOPTIONS` and `ACQUIPARAMS` structures described in following chapters.


A high speed video camera temporarily stores recorded movies in the internal video memory, called *Frame Buffer Memory (FBM)*. The frame buffer memory can be split into multiple partitions to quickly record successive events in different movies. The memory of each partition is managed to work like a circular buffer: during a long recording, even if memory gets filled, the acquisition will continue by overwriting the oldest images from that partition.

Each partition has three main buffers allocated for acquisition parameters, image processing parameters and other parameters. The first one allows configuring recording parameters. As soon as the acquisition ends, these parameters become read-only.


Image processing parameters are not affecting recorded image data. They are stored as meta-data and used when images are displayed on video output or on PC applications. Unlike the acquisition parameters, they cannot be preconfigured. They are inherited from one special camera structure when the recording of the cine ends. This camera structure is holding a set of global image processing parameters used for processing the images of the movie being currently recorded. As soon as the acquisition ends, values from the camera data structure are copied into the image processing buffer of the partition the movie was recorded on. From this point forward, the image processing parameters of the

recorded movie can be changed as needed by directly accessing its own image processing buffer.

Besides the main partitions where recorded movies are stored, there is a hidden video memory partition used for preview purposes. If all partitions contain finished recordings, this is where the camera will store one live image. In this case, the camera uses a smaller acquisition frame rate thus reducing power consumption. This working mode is called *pretrigger mode*. A camera automatically switches to *pretrigger mode* when all main partitions become full. A camera can also be set in *pretrigger mode* even if the FBM is empty; in this mode one can preview images for a certain set of acquisition parameters and, in the same time, reduce power consumption.

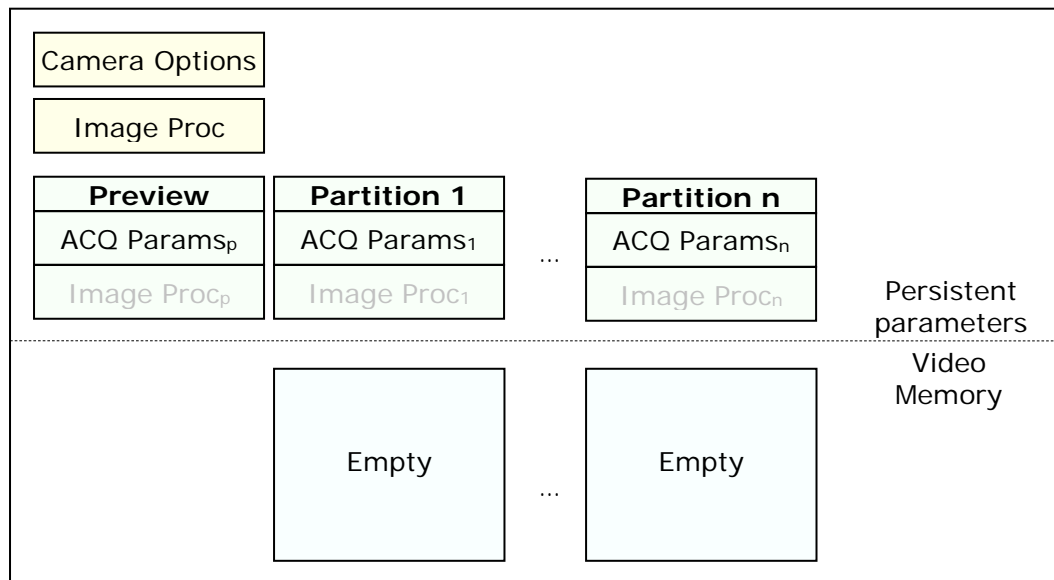

**Figure 2-1. Camera data organization before starting a new acquisition**

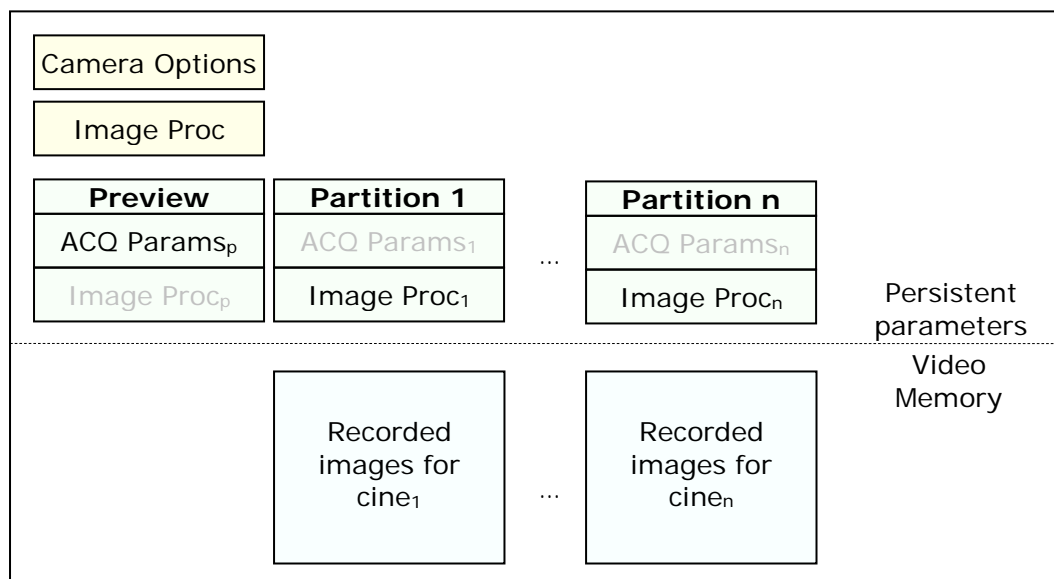

**Figure 2-2. Camera data organization after finishing the acquisition**

### 2.1.2. The Cine

A *cine* is a generalization of a Phantom camera movie. It is a collection of one or more stored images accompanied by a header describing the cine. The header contains information such as image resolution, acquisition frame rate or image processing parameters. A *cine* allows a unitary, simplified approach no matter if working with a movie from the latest generation Vision Research camera or from older camera models. It also enables similar handling over different types of movie storage (the movie can reside on the camera, on special storage devices, like the CineMag, or it can be stored in a file).

While most of the operations are performed in the same manner for all types of cines, a cine storage type imposes some distinctions over what actions are allowed or not. For example, only a recorded cine stored inside the RAM of the camera can be saved to a CineMag; a file stored cine movie cannot be saved to a CineMag. Thus, the two main types of cines are: camera *cines* and file *cines*.

A camera *cine* refers to a cine for which images and most of the header data are stored on camera. Such a cine can be a recorded movie stored either in camera's RAM or in special storage devices like the CineMag.

A cine stored in the RAM of the camera is volatile; it will be lost if the camera is turned off or unpowered. A cine stored in the camera flash (like the CineMag) is nonvolatile; it will not be lost if the camera is turned off.

A file *cine* generally refers to a cine for which images and header data are stored in file. It can be one of Vision Research proprietary formats (*.cine*, *.cci*) or any other video formats the SDK offers reading support for. Check Supported File Formats for more details.

Phantom SDK offers a way to perform basic cine operations in the same way no matter the cine type, by introducing another special element: the *CINEHANDLE*. A *CINEHANDLE* is a handle to a cine. It facilitates access to cine header or image data and it enables simple control for actions such as saving or converting. Most of the cine operations allowed by the Phantom SDK require a *CINEHANDLE* object to identify the cine.

Remember to always use a *CINEHANDLE* object for reading images from a cine!

Please note that the two notions, the cine – as an abstract generalization of a Vision Research high speed camera movie – and the *CINEHANDLE* – as a programming element used for referring to a cine – may overlap in this document. For example, a *CINEHANDLE* input parameter of a function may be named the input cine.

### 2.1.3. Image Numbering

The images are numbered in a growing order using 32 bit signed values. The images before trigger have negative numbers. The images after the trigger have zero or positive numbers.

## 2.2. Relations

### 2.2.1. The Camera Pool

All Phantom cameras the Phantom SDK established connections with at a given time are referred to as the *camera pool*. If a camera is not in the *camera pool* than the SDK is not aware of its presence, cannot communicate with it, cannot read images or set any of its acquisition parameters.

The *camera pool* is actually a buffer to keep information about connected cameras. Inside the camera pool each camera is identified based on a numerical index, called the *camera number*. The first index is 0. The *camera number* is used a lot in Phantom SDK when controlling a camera or reading/writing its parameters.

Building a camera pool can be done by setting the software to discover all cameras that are available on the network or by entering a list of the IP addresses of the cameras we want to connect to.

In the first case a "discovery protocol" is used. The computer sends a broadcast message on the network and the cameras answer to it providing their IP address.

An important feature for a camera control system is to be able to discover new added cameras and to support the disconnecting and reconnecting of the cameras.

Before version 770, applications based on Phantom SDK did that by discovering changes in camera pool (calling `PhCheckCameraPool`) and updating camera pool (calling `PhNotifyDeviceChangeCB`). Last call can be done only when there are no pending operation on cameras.

The enumeration is now fully controlled inside the SDK. Application call `PhConfigPoolUpdate` to enable the automated refresh of the camera pool.

New cameras that appear during the application run could be added at the end of the camera pool, so the existing cameras will not change number and the communication to them will not be affected at all.

The application will notice a new camera is available just by the increase of the number returned by `PhGetCameraCount`.

If the communication to an existing camera is broken (a timeout occurred) the software will see camera in a special "Offline" status. If the camera gets back online later, the connection will be restored and camera will be "revived" in the same position in the pool.

Existing cameras that go Offline will be found by the Timeout result at the commands sent to the camera. Both new cameras on network and the Offline revived cameras will be found by the network discovery protocol based on broadcast and UDP responses.

### 2.2.2. Cameras and Cines

As their name implies, camera cines belong to a camera. If a camera is identified based on a camera number, a camera cine is identified based on:

- the number of the camera it belongs to
- a *cine number*

The cine number represents the partition number the cine is stored on or recorded on. The first partition number is 1. For example, if a camera has three partitions, corresponding cines numbers will be: 1, 2, and 3.

Special values are used for cines stored in cameras NVMEM (flash or CineMag). They are indexed beginning with the value returned by the function `PhFirstFlashCine`.

While a camera cine as previously defined is strictly bound to one of the partitions of a camera, there is a special type of camera cine that does not seem to follow the rules. However, it is extremely useful and necessary. This camera cine is called *live cine*. The term *current cine* may also be used.

The *live cine* is the cine the camera is currently recording images into, thus it is not necessarily associated with the same partition number at different moments in time. Each camera has exactly one *live cine*.

An application will often need to display images as currently recorded by the camera. When requesting such images, one cannot precisely tell from which partition number or at which resolution these will be taken. This happens because from the moment the image request is sent and the moment the camera processes it, the current cine number and acquisition parameters may change. Here is where the *live cine* becomes most useful: an application can use a *live cine* entity for reading images like any other cine letting the SDK and the camera to solve any details about the actual cine the camera is recording images into.

As previously mentioned, for basic cine access and control one should use a *CINEHANDLE* object. The same stands true for camera cines no matter if they are stored cines or live cines. While the camera and cine number physically identify a camera cine and also allow low level control, a *CINEHANDLE* object will greatly simplify operations such as reading images, reading time information, setting image processing parameters or saving to file. Except for the creation step of a *CINEHANDLE* object, all of the above actions are performed in the same way for either a live cine, a camera recorded cine or a file stored cine.

Please note that while a cine refers to a set of recorded images plus a header, the term *cine number* may sometimes be used in this documentation to identify the partition where a new cine will be recorded.

## 2.3. Actions

### 2.3.1. Camera and Cine Control

Phantom SDK functions also allow controlling camera and cine objects. Most camera operations directly refer to camera capabilities and features and are not discussed in detail in this document. Typical camera actions include but are not limited to:

- operations for preparing a new recording session
- recording control operations
- playing back recorded movies
- saving recorded movies to the computer hard disk
- saving recorded movies to camera flash memory (like the CineMag)
- configuring video output parameters
- playing back cine images on video output or recording this output to tape devices

Main cine control operations are related to cine image access, image processing and cine saving operations. A camera cine stored either in camera RAM or on



special storage devices (like the CineMag) can be saved to a file using our proprietary cine file format or any other supported file types (see Supported File Formats chapter). Phantom SDK allows image related operations for file cines as well. Similar to a camera cine, a file cine can be saved to various cine file formats. Saving a file cine to another file is also referred to by the term *conversion*.

### 2.3.2. Cine Use Case

As previously defined, a cine is a collection of images plus a header describing the images and also containing different metadata. An important part of this header is represented by the image processing parameters. Their values are stored separately so that the original sensor data (raw pixels) remains untouched allowing to "undo" any processing.

An application may need to display images from a cine, perform different analysis and measurement operations, save one or more images of this cine to various file formats. For each of these actions, images of the input cine are read, processed and further displayed or written to a file. However, none of these operations will affect the input raw cine. The input raw cine remains a collection of images accompanied by a header. For example, requesting images from a color camera CINE in order to be displayed will involve a demosaicing process and will produce images with 3 samples per pixel. Nevertheless, the source cine remains a cine containing images with 1 sample per pixel. The cine is not affected by any operations performed on its images after being read.

In order to indicate more intuitively all processing that must be done for images of the input cine, the user will specify a *Use Case* for that cine. It is a simple way of saying what a user wants to do with a cine, with its images. One may want to get cine images in order to later display them or to save them in one of the supported file formats.

It is easier to specify that a certain set of camera cine images are to be displayed (hence color interpolated) or to be saved to a raw format (no color interpolation) than to explicitly call each processing one by one.

Phantom SDK introduced use cases as a way of simplifying image processing related configuration. Together with options that allow customizing use cases, it offers a very flexible and easy to use solution for either basic or advanced cine image processing.

### 2.3.3. Putting It All Together - Creating a New Application

Any application built upon Phantom SDK should follow these basic steps:

- Register the application as a Phantom SDK client
  - ✓ `PhRegisterClientEx`
- Obtain a `CINEHANDLE` for the live cine
  - ✓ `PhGetCineLive`
- Prepare for a new recording session by configuring the camera and setting-up acquisition parameters for each partition
  - ✓ `PhSetCameraOption`
  - ✓ `PhSetSingleCineParams/PhSetCineParams`
  - ✓ `PhSet`
- In the same time, display live images
  - ✓ `PhGetCineImage`
- Start the recording
  - ✓ `PhRecordCine`



- End the recording
  - ✓ PhSendSoftwareTrigger
- Obtain a CINEHANDLE for the newly registered cine
  - ✓ PhNewCineFromCamera
- Display cine images
  - ✓ PhGetCineImage
- Adjust recorded cine image parameters
  - ✓ PhSetCineInfo
- Read recorded cine acquisition parameters
  - ✓ PhGetCineParams
  - ✓ PhGetCineInfo
  - ✓ PhGet
- Play images from the recorded cine
  - ✓ PhVideoPlay
- Finalize
  - ✓ PhDestroyCine
- Unregister
  - ✓ PhUnregisterClientEx

### 3. Data Structure Reference

#### 3.1. Time Related Data Structures

##### 3.1.1. TIME64

```
typedef UINT32 FRACTIONS, *PFRACTIONS;

typedef struct tagTIME64
{
    FRACTIONS fractions;
    UINT32 seconds;
} TIME64, *PTIME64;
```

Header: PhInt.h, PhCon.h

##### *fractions*

Fractions of a second.

Stored here multiplied by  $2^{32}$  and rounded to integer.

Least significant 2 bits store information about IRIG synchronization and the camera Event input:

Bit0 Value	Meaning
0	IRIG synchronized
1	Not synchronized

Bit1 Value	Meaning
0	Event input = 0 (short to ground)
1	Event input = 1 (open)

##### *seconds*

Seconds from January 1, 1970.

Compatible with the C library routines.

The maximum year allowed by the unsigned representation is 2106.

**TIME64** is a compact format for storing very precise and absolute time. The internal representation is fixed point 32.32, the unit is second and the origin is January 1, 1970. The resolution of this representation is below 1 nanosecond and the allowed interval of years is from 1970 to 2106.

The least significant bit of *fractions* stores the information whether the IRIG is locked to an external signal. The second least significant bit (bit 1) stores the status of Event, a general-purpose binary input.

The type of the *seconds* field was changed to `UINT32` (from `time_t`) to enhance the compatibility with Visual Studio 2005 and above while maintaining its size to 32 bits.

### 3.1.2. TC

```
typedef struct tagTC
{
    unsigned char framesU:4;
    unsigned char framesT:2;
    unsigned char dropFrameFlag:1;
    unsigned char colorFrameFlag:1;
    unsigned char secondsU:4;
    unsigned char secondsT:3;
    unsigned char flag1:1;
    unsigned char minutesU:4;
    unsigned char minutesT:3;
    unsigned char flag2:1;
    unsigned char hoursU:4;
    unsigned char hoursT:2;
    unsigned char flag3:1;
    unsigned char flag4:1;
    UINT userBitData;
}TC, *PTC;
```

Header: PhInt.h, PhCon.h

#### *framesU*

Units of frames.  
4 bits are used for storing this value.

#### *framesT*

Tens of frames.  
2 bits are used for storing this value.

#### *dropFrameFlag*

Dropframe flag.  
Currently not used. It is always 0.

#### *colorFrameFlag*

Colorframe flag.  
Currently not used. It is always 0.

#### *secondsU*

Units of seconds.  
4 bits are used for storing this value.

#### *secondsT*

Tens of seconds.  
3 bits are used for storing this value.

#### *flag1*

Currently not used. It is always 0.

#### *minutesU*

Units of minutes.  
4 bits are used for storing this value.

#### *minutesT*

Tens of minutes.  
3 bits are used for storing this value.

#### *flag2*

Currently not used. It is always 0.

#### *hoursU*

Units of hours.  
4 bits are used for storing this value.

#### *hoursT*

Tens of hours.

2 bits are used for storing this value.

*flag3*

Currently not used. It is always 0.

*flag4*

Currently not used. It is always 0.

*userBitData*

32 bits for storing user data.

Currently not used.

**TC** is used for representing time code according to SMPTE 12M-1999 standard. Time data, flag bits and user bits are all represented in **TC** precisely following the standard data packing specifications.

A time code is a 64bit code assigned to each video frame consisting of hours, minutes, seconds, a frame number, all counting for 32 bits, plus another 32 bits for storing user data. The time code is usually limited to one day duration. The frame number goes from 0 to (playback rate – 1).

The hours, minutes, seconds, frames are represented using BCD. From the four bits reserved for each digit some may never be used. The tens of hours, for example, cannot be larger than 2. Such unnecessary bits are assigned special uses (flags).

Note that the maximum playback rate is 40 (the maximum tens of frames is 3).

### 3.1.3. TCU

```
typedef struct tagTCU
{
    UINT frames;
    UINT seconds;
    UINT minutes;
    UINT hours;
    BOOL dropFrameFlag;
    BOOL colorFrameFlag;
    BOOL flag1;
    BOOL flag2;
    BOOL flag3;
    BOOL flag4;
    UINT userBitData;
}TCU, *PTCU;
```

**Header:** PhInt.h, PhCon.h

*frames*

Frame index.

*seconds*

Number of seconds.

*minutes*

Number of minutes.

*hours*

Number of hours.

*dropFrameFlag*

Dropframe flag.

Currently not used. It is always 0.

*colorFrameFlag*

Colorframe flag.

Currently not used. It is always 0.

*flag1*

Currently not used. It is always 0.

*flag2*

Currently not used. It is always 0.  
*flag3*  
 Currently not used. It is always 0.  
*flag4*  
 Currently not used. It is always 0.  
*userBitData*  
 32 bits for storing user data.  
 Currently not used.

**TCU** is used for representing SMPTE time code. However, unlike **TC**, it does not follow the standard specifications regarding data packing. The *frames* index, the *hours* value, the *minutes* value, the *seconds* value are each represented by a decimal coded value.

## 3.2. Acquisition Parameters

### 3.2.1. ACQUIPARAMS

```
typedef struct tagACQUIPARAMS
{
    UINT ImWidth;
    UINT ImHeight;
    UINT FrameRate;
    UINT Exposure;
    UINT EDRExposure;
    UINT ImDelay;
    UINT PTFrames;
    UINT ImCount;
    UINT SyncImaging;
    UINT AutoExposure;
    UINT AutoExpLevel;
    UINT AutoExpSpeed;
    RECT AutoExpRect;
    BOOL Recorded;
    TIME64 TriggerTime;
    UINT RecordedCount;
    INT FirstIm;
    UINT FRPSteps;
    INT FRPImgNr[16];
    UINT FRPRate[16];
    UINT FRPExp[16];
    UINT Decimation;
    UINT BitDepth;
    UINT CamGainRed;
    UINT CamGainGreen;
    UINT CamGainBlue;
    UINT CamGain;
    BOOL ShutterOff;
    UINT CFA;
    char CineName[256];
    char Description[4096];
    UINT FRPShape[16];
    double dFrameRate
} ACQUIPARAMS, *PACQUIPARAMS;
```

Header: PhCon.h

#### Image Format

*ImWidth*

Image width.

*ImHeight*

Image height.

*BitDepth*

Acquisition bit depth.

Read-only field.

*CFA*

Color Filter Array of the camera the cine was recorded on.

*CamGainRed*

Red gain attached to a cine saved in the magazine. Normally it tells the White balance at recording time.

*CamGainGreen*

Green gain attached to a cine saved in the magazine.

*CamGainBlue*

Blue gain attached to a cine saved in the magazine.

*CamGain*

Global gain.

## Main Acquisition Parameters

*FrameRate*

Frame rate in frames per second.

It is the requested frame rate represented as an integer. The camera will realize the closest possible frame rate, depending on its internal clock frequency. This parameter is normally constant during recording.

However, it is possible to dynamically change the frame rate by using the frame rate profile. See *FRPSteps*, *FRPImgNr*, *FRPRate* fields.

The actual frame rate can be obtained from the difference of two successive image times. To read a cine image time use *PhGetCineAuxData* with the *GAD\_TIME* selector.

*dFrameRate*

High precision acquisition frame rate, replace *uint32\_t FrameRate*

*Exposure*

Exposure duration in nanoseconds.

It is the requested value; the camera will do the closest possible value.

If the autoexposure mode is enabled, the exposure can change dynamically during the recording. See *AutoExposure*, *AutoExpLevel*, *AutoExpRect* fields.

To get the actual exposure duration for any acquired cine image use *PhGetCineAuxData* with the *GAD\_EXPOSURE* selector.

*PTFrames*

Count of frames to be recorded after trigger.

It can be any non-zero, positive integer, up to 4 billion.

A value greater than *ImCount* has the effect of an adjustable delay; the frame buffer is a circular buffer and thus only the last *ImCount* images will be stored here.

*SyncImaging*

Defines the source of image acquisition pulses. Possible values are:

Value	Meaning
<i>SYNC_INTERNAL</i>	Repetitive acquisition at the specified <i>FrameRate</i> frequency
<i>SYNC_EXTERNAL</i>	It starts the exposure of each frame when an external signal rises
<i>SYNC_LOCKTOIRIG</i>	It allows the synchronization of multiple cameras using the common IRIG time signal

## Cine Status Information

### *Recorded*

If this field is `TRUE`, the cine was recorded and is available for playback. *TriggerTime*, *RecordedCount*, *ImCount*, *FirstIm* fields contain valid and final values. They are now read only fields.

### *CineName*

Cine name.

### *Description*

Cine description on maximum 4095 chars.

### *TriggerTime*

Exact trigger time for a recorded cine. Its value is somewhere between the time of image -1 and the time of image 0.

Read-only field.

### *ImCount*

For an empty partition, this field represents the maximum number of images that can be stored using current acquisition parameters.

For a partition holding a recorded cine, this field represents the number of cine images actually stored on this partition.

Read only field.

### *RecordedCount*

Number of recorded frames.

*ImCount* may be smaller if the cine was partially saved to flash memory.

Read-only field.

### *FirstIm*

First image number of this cine; it may be different from *PTFrames - ImCount* if the cine was partially saved to Non-Volatile memory.

Read-only field.

## Advanced Acquisition Parameters

### *EDRExposure*

Extreme dynamic range exposure duration (nanoseconds).

### *ImDelay*

Delay in nanoseconds between the sync/irig and the start of the image exposure.

It is not available in all camera models.

### *AutoExposure*

Control of the exposure duration based on the light conditions.

Value	Meaning
0	Disabled
1	Enabled, lock at trigger
3	Enabled, active after trigger

If auto exposure mode is enabled, the camera continuously adjusts the exposure duration to get the best images, even when light conditions change. In this case, the exposure the camera chooses is never greater than the *Exposure* value.

### *AutoExpLevel*

Autoexposure level to be realized by the automatic exposure system.

Range of possible values is [0...255].

### *AutoExpSpeed*

Speed for autoexposure control.

Currently not used.

### *AutoExpRect*

Rectangle representing the region of interest where the autoexposure measures the average pixel level.

### *ShutterOff*

This BOOL setting puts the shutter off to force the maximum exposure of the camera.

It is useful for PIV where the camera is expected to have the shortest interval between successive exposures.

## Frame Rate Profile

### FRPSteps

Number of frame rate profile steps:

Value	Meaning
0	No frame rate profile
any other value	Number of points where the frame rate changes

The frame rate can be changed only after the trigger moment and the change can be a step to another value that will remain constant (flat) or can be continuously variable from a value to another value (ramp).

### FRPImgNr

Image number where to change the acquisition frame rate.

The frame rate can be changed only after the trigger moment.

The acquisition parameters used before *FRPImgNr[0]* are the pretrigger parameters as stored in *FrameRate* and *Exposure* fields.

Allocated for 16 points.

### FRPRate

New value for frame rate (fps).

At image number *FRPImgNr[i]* the frame rate becomes *FRPRate[i]*.

Allocated for 16 points.

### FRPShape


Frame rate evolution shape during recording.

Read-only parameter.


At image number *FRPImgNr[i]* the frame rate becomes *FRPRate[i]* following a shape described by *FRPShape[i]*.

All the steps share the same shape: ramp in cameras that support it and flat in the older cameras.

Value	Meaning
0	flat
1	ramp

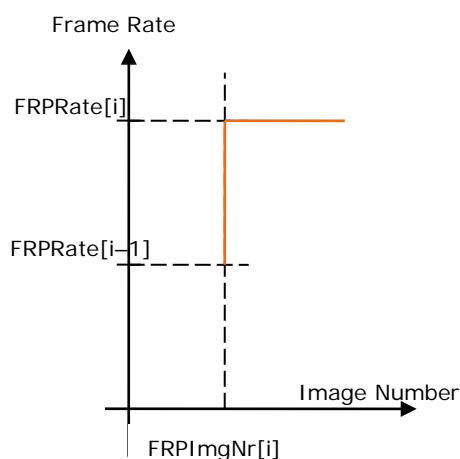

Flat Shape

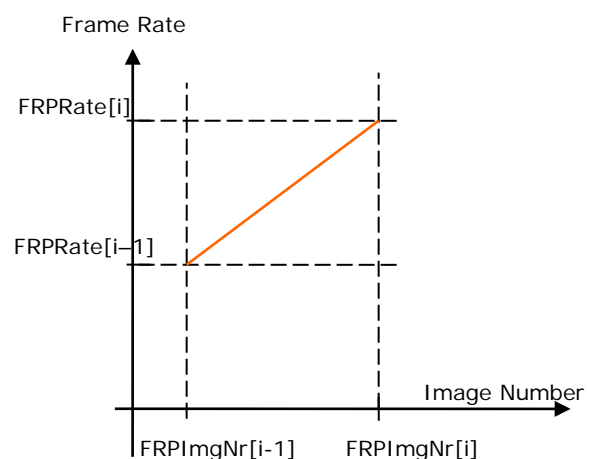

Ramp Shape

### FRPExp

New value for exposure (nanoseconds).



Currently not implemented in cameras - the exposure cannot be changed during the profile. All frame rates specified in the profile have to be smaller than the reciprocal of the exposure.

## Other

### Decimation

Factor by which to reduce the frame rate when sending images from camera. Currently used for i3 external memories.

The **ACQUIPARAMS** structure gathers most of the camera parameters related to the acquisition process from basic acquisition parameters such as exposure time to more advanced options like frame rate profile. It is used when configuring a camera for a new acquisition, as well as when reading the parameters, a cine was recorded with.

To configure a camera for a new recording, use `PhSetCineParams` or `PhSetSingleCineParams` functions depending on whether the camera is set to work on multi cine or single cine mode. To read the acquisition parameters from a stored camera cine use `PhGetCineParams` function.

The acquisition parameters as depicted by the **ACQUIPARAMS** structure are also stored in the nonvolatile memory of the camera. This way, the last settings remain active even after the camera is powered off and restarted.

Newer acquisition parameters are controlled by calling `PhSet/PhGet` functions with a corresponding selector (see Appendix 3). These are also the functions future Phantom SDK versions will develop for acquisition parameters control.

## 3.3. Cine Related Data Structures

### 3.3.1. IMRANGE

```
typedef struct tagIMRANGE
{
    INT First;
    UINT Cnt;
} IMRANGE, *PIMRANGE;
```

Header: PhCon.h

#### *First*

First image number.

#### *Cnt*

Count of images.

The **IMRANGE** defines a range of images by specifying the first image and the image count.

### 3.3.2. CINESTATUS

```
typedef struct tagCINESTATUS
{
    BOOL Stored;
    BOOL Active;
    BOOL Triggered;
} CINESTATUS, *PCINESTATUS;
```

Header: PhCon.h

#### *Stored*

Recording of this cine finished.

*Active*

Acquisition is currently taking place in this cine.

*Triggered*

The cine has received a trigger.

The cine status structure provides information about whether a recording finished or the trigger signal was received. Use `PhGetCineStatus` to obtain a **CINESTATUS** array with status information for all the cines in a certain camera.

## 3.4. Camera Related Data Structures

### 3.4.1. CAMERAOPTIONS

```
typedef struct tagCAMERAOPTIONS
{
    UINT OSD;
    UINT VideoSystem;
    BOOL ModulatedIRIG;
    BOOL RisingEdge;
    UINT FilterTime;
    BOOL MemGate;
    BOOL StartInRec;
    BOOL Color;
    BOOL TestImage;
    RGBQUAD OSDColor[8];
    BOOL OSDOpaque;
    int OSDLeft;
    int OSDTop;
    int OSDBottom;
    int ImageX;
    int ImageY;
    BOOL AutoSaveNVM;
    BOOL AutoSaveFile;
    BOOL AutoPlay;
    BOOL AutoCapture;
    CHAR FileName[MAX_PATH];
    IMRANGE SaveRng;
    BOOL LongReady;
    BOOL RealTimeOutput;
    UINT RangeData;
    UINT SourceCamSer;
    UINT SliceNr;
    UINT SliceCnt;
    BOOL FRPi3Trig;
    UINT UT;
    int AutoSaveFormat;
    UINT SourceCamVer;
    UINT RAMBitDepth;
    int VideoTone;
    int VideoZoom;
    int FormatWidth;
    int FormatHeight;
    UINT AutoPlayCnt;
    UINT OSDDisable;
    BOOL RecToMag;
    BOOL IrigOut;
    int FormatXOffset;
    int FormatYOffset;
} CAMERAOPTIONS, *PCAMERAOPTIONS;
```

**Header:** PhCon.h

## On Screen Display Options

### *OSD*

If TRUE, it enables the camera On Screen Display.

### *Color*

If TRUE, it enables OSD color in PAL/NTSC modes.

### *OSDColor*

Array of 8 RGBQUAD elements describing various OSD colors. They have the following meaning:

Array Element	Meaning
<i>OSDColor[0]</i>	Background color for OSD fields.
<i>OSDColor[1]</i>	Logo color in waiting for trigger mode.
<i>OSDColor[2]</i>	Logo color in trigger mode.
<i>OSDColor[3]</i>	Logo color in cine stored mode.
<i>OSDColor[4]</i>	Logo color in preview mode.
<i>OSDColor[5]</i>	Foreground color for camera identification fields.
<i>OSDColor[6]</i>	Foreground color for acquisition parameter fields.
<i>OSDColor[7]</i>	Foreground color for status fields.

This parameter may not be available in newer cameras.

### *OSDOpaque*

OSD text is opaque or transparent.

### *OSDLeft*

Extra border on the left of the OSD text.

### *OSDTop*

Extra border above the OSD text.

### *OSDBottom*

Extra border below the OSD text.

### *UT*

Time format as displayed by the OSD:

Value	Meaning
0	Local time.
1	Universal Time (GMT).

### *OSDDisable*

Selective disable of the analog and digital OSD:

Value	Meaning
0	OSD on.
1	OSD turned off on monitor.
2	OSD turned off on recorder output.
3	OSD turned off on both the monitor and the recorded output.

## Production Area Rectangle

### *FormatWidth*

Width for the rectangle to overlap on the video image.

### *FormatHeight*

Height for the rectangle to overlap on the video image.

### *FormatXOffset*

x offset for the rectangle to overlap on the video image.

### *FormatYOffset*

y offset for the rectangle to overlap on the video image.

## Video Out Control

### *VideoSystem*

Video output standard used.

Possible values are:

0	NTSC
1	PAL
4	720P60
5	720P50
12	720P59
20	1080P30
21	1080P25
28	1080P29
36	1080P24
44	1080P23
52	1080PSF30
53	1080PSF25
60	1080PSF29
68	1080I30
69	1080I25
76	1080I29
84	1080PSF24
92	1080PSF23

### *TestImage*

If TRUE, a color bar test image is displayed on the connected monitor.

### *ImageX*

Image X coordinate on screen.

### *ImageY*

Image Y coordinate on screen.

### *AutoPlayCnt*

Repeat count for the auto playback to video.

Value	Meaning
0	If <i>AutoCapture</i> is 0, the video play back loops forever. If <i>AutoCapture</i> is non zero, the video is played back only once and afterwards the acquisition is restarted.
Non zero value	Loop counter for the video play back.

### *VideoZoom*

Zoom of the video image:

Value	Meaning
0	Fit.
1	Zoom 1.

## Global Camera Options

### *RisingEdge*

Trigger polarity.

Value	Meaning
TRUE	Camera triggers on the rising edge of the trigger signal.
FALSE	Camera triggers on the falling edge of the trigger signal.

### *ModulatedIRIG*

If TRUE, IRIG input accepts modulated signal.

### *FilterTime*

A time constant related to the filter the camera uses for ignoring spurious trigger signals.

Range of possible values: [ 0 . . 255 ].

### *MemGate*

Specifies the function for one of the capture connector pins.

For this pin, a TTL supplied input signal can be used to toggle the Memory Gate or as a Pre-trigger input signal.

When the Memory Gate feature is enabled, an input signal, when activated, prevents the storing of the acquired image data by disabling the write-access to memory.

A pre-trigger input, on the other hand, allows the user to decide when to place the camera into the Ready-State or Capture Mode.

Value	Meaning
TRUE	Memgate function
FALSE	Pre-trigger function

#### *StartInRec*

Specifies the default operational state a camera should be placed into when starting up.

Value	Meaning
TRUE	Start in recording mode. The camera is recording images and is waiting for trigger.
FALSE	Start in preview mode. The camera is waiting for pretrigger to start the recording. This mode saves power.

#### *IrigOut*

If TRUE, change the Strobe pin to Irig Out at Miro cameras.

#### *LongReady*

Value	Meaning
TRUE	Ready signal is 1 from the start to the end of recording.
FALSE	Ready signal is 1 from the start of recording until the cine is triggered.

#### *RangeData*

Range data size per image in bytes.

Value	Meaning
0	No range data.
Any other value	The range data size in bytes per image.

#### *RecToMag*

If TRUE, acquired images are recorded directly into the Cine Mag.

#### *RealTimeOutput*

If TRUE, it enables the real time output.

### **Camera Autosave**

#### *AutoSaveNVM*

If TRUE, a camera having an internal flash storage will save the just-recorded cine to this storage.

#### *AutoSaveFile*

If TRUE, the camera will save the just-recorded cine to a file specified by the *FileName* field.

#### *AutoPlay*

If TRUE, the camera will play back the recorded cine on video immediately after the recording ends and after the autosave operations, if any, are finished.

#### *AutoCapture*

If TRUE, the camera will restart the acquisition immediately after the current recording ends and after the autosave operations, if any, are finished.

#### *FileName*

File path where to save a cine if *AutoSaveFile* is enabled.

### *SaveRng*

Image range for the above operations.

### *AutoSaveFormat*

Frame format used for file/flash cine save operations:

Value	Meaning
8	Calibrated 8bit format.
16	Calibrated 16bit format.
266	Calibrated 10bit packed format. Only for CineMag.

## Image Format

### *RAMBitDepth*

Bit depth used by the camera when acquiring image data.

### *VideoTone*

It selects a predefined tone curve for video.

## Image Cube

### *SourceCamSer*

Source camera serial.

### *SliceNr*

Image Cube Slice number (0, 1, 2 ...).

### *SliceCnt*

Total count of Image Cube slices connected to a certain camera.

### *FRPi3Trig*

Frame rate profile start at Image Cube trigger.

### *SourceCamVer*

Source camera version (external memory slice).

The **CAMERAOPTIONS** structure is used for general camera configuration: video output control, parameters adjustment for the OSD, camera behavior when restarting or when an acquisition ends.

To read current camera configuration, use `PhGetCameraOptions` function.

To change any of the **CAMERAOPTIONS** fields, read the current settings, update the fields you want to change and write the updated structure back to the camera by calling `PhSetCameraOptions` function.

### 3.4.2. CAMERAID

```
typedef struct tagCAMERAID
{
    UINT IP;
    UINT Serial;
    char Name[256];
    char Model[256];
} CAMERAID, *PCAMERAID;
```

Header: PhCon.h

#### *IP*

Integer value representing a camera IP address.

Each of the four IP address values is stored on one unsigned byte.

#### *Serial*

Camera serial number.

#### *Name*

Camera name.

#### *Model*

Camera model.

**CAMERAID** is a simple structure storing basic information about a camera. It is currently used by `PhGetAllIpCameras` function.

### 3.4.3. PULSEPARAM

```
typedef struct tagPULSEPARAM
{
    UINT Size;
    BOOL Invert;
    BOOL Falling;
    double Delay;
    double Width;
    double Filter;
} PULSEPARAM, *P PULSEPARAM;
```

**Header:** PhCon.h

#### *Size*

Size of the structure, for versioning in the future if needed

#### *Invert*

Invert the signal at the output of the pulse processor

#### *Falling*

Selects falling edge mode for the pulse processor. This mode is only relevant if width is also specified. When the falling token is present together with width, the pulse processor will generate a negative pulse, triggered from the negative edge of the input signal

#### *Delay*

Delay the output pulse by the specified time. clocks. If the width token is not present, both edges of the signal are delayed by the same amount. If the width token is present, the delay is measured from the rising edge of the input signal unless the falling token is present, in which case the delay is measured from the falling edge of the input. The delay time is specified in seconds and is internally converted and rounded to pixel clock units. The maximum delay time is at least 10 seconds

#### *Width*

When a width token is present, a defined-length pulse is generated, which starts after the specified delay after the active edge of the input signal. The length of the pulse is specified in seconds and internally converted to pixel clock units. The maximum pulse width is at least 10 seconds. However, if the period of the input signal is lower than the width, the latter is dynamically clamped to the signal period. The minimum pulse width is one pixel clock

#### *Filter*

When a filter token is present, the input of the pulse processor is filtered through an edge filter of the specified time. In order for the output of the filter to be asserted, the input signal must be continuously asserted for the specified duration. In order for the output of the filter to be deasserted, the input signal must be continuously deasserted for the same duration. The edges of the input are thus delayed by the specified filter time (for a "clean" input pulse). Filtering is applied before and independently of the delay and duration. The filter time is specified in seconds and is internally converted and rounded to pixel clock units. The maximum filter time is 1 second

**PULSEPARAM** is a simple structure for control of programmable I/O operation. It is currently used by PhGet1/PhSet1 function.

## 3.5. Save Related Data Structures

### 3.5.1. SAVECINEINFO

```
typedef struct tagSAVECINEINFO
{
    CINEHANDLE CH;
    UINT progressPercent;
}
SAVECINEINFO, *PSAVECINEINFO;
```

Header: PhFile.h

*CH*

Cine handle.

*progressPercent*

Current progress for a write cine to file operation.

**SAVECINEINFO** is used to store information about a write cine file in progress operation.

## 3.6. Image Headers

Two data structures are used by the Phantom SDK to describe an image:

**BITMAPINFOHEADER** and **IH**.

**BITMAPINFOHEADER** is the Windows data structure containing size and color image information. However, the Phantom SDK extends some fields meaning in order for images with more than 8bits per color component to be supported.

**IH** is the Phantom SDK specific image header, built as an extension of the **BITMAPINFOHEADER** structure.

Even if a **IH** structure is an extension of a **BITMAPINFOHEADER** structure, be careful to never use a cast for converting from one type to the other. Where one or more images are described by a **IH** structure while a **BITMAPINFOHEADER** one is needed, call **PhIHtoBITMAPINFOHEADER** function.

### 3.6.1. BITMAPINFOHEADER

```
typedef struct tagBITMAPINFOHEADER
{
    DWORD    biSize;
    LONG     biWidth;
    LONG     biHeight;
    WORD     biPlanes;
    WORD     biBitCount;
    DWORD    biCompression;
    DWORD    biSizeImage;
    LONG     biXPelsPerMeter;
    LONG     biYPelsPerMeter;
    DWORD    biClrUsed;
    DWORD    biClrImportant;
} BITMAPINFOHEADER, *PBITMAPINFOHEADER;
```

Header: PhInt.h



*biSize*

The number of bytes required by this structure.

*biWidth*

Image width in pixels.

*biHeight*

Image height in pixels.

Phantom SDK specific:

This field is always positive. For BI\_RGB formats, images are always stored bottom up. For BI\_PACKED formats, images are always stored top-down.

*biPlanes*

It specifies the number of planes for the target device.

This value must be set to 1.

*biBitCount*

It specifies the number of bits-per-pixel used for storing a pixel value.

Phantom SDK specific:

Possible *biBitCount* values are limited to 8 and 16 for monochrome images, 24 and 48 for RGB color images.

The meaning of a 16 bit image is different from Windows: it is a 16 bit per pixel gray image. Each pixel value is stored on 16 bits even if the real bit depth produced by the camera is 14, 12 or 10 bits. 48 value of this field corresponds to a color image having 16 bits per color component.

Color palette images (8 bpp) are not supported; in the Phantom environment they are converted to 24 bpp after the file load or after the copy from clipboard. The palette is not written in the cine file but a gray palette is needed to render the monochrome 8bpp DIBs in Windows.

*biCompression*

It specifies the type of compression for a compressed bottom-up bitmap (top-down DIBs cannot be compressed).

Phantom specific:

Only BI\_RGB = 0 and BI\_PACKED = 256 are supported.

When BI\_PACKED is used, image data is packed without using any padding bits. For example, for 10bit image data, 4 sample values are stored in 5 bytes. This packing reduces used memory or file size by a factor of 10/16. A BI\_PACKED bitmap is stored in top-down row order both in a packed cine file and in the camera memory.

*biSizeImage*

Image size in bytes.

*biXPelsPerMeter*

It specifies the horizontal resolution, in pixels-per-meter, of the target device for the bitmap. An application can use this value to select a bitmap from a resource group that best matches the characteristics of the current device.

Phantom specific:

*biXPelsPerMeter*, *biYPelsPerMeter* are resolutions computed at the camera sensor level. To get the scene resolution, you have to multiply these values by the distance from the camera to the scene divided by the focal length of the lenses.

*biYPelsPerMeter*

Vertical resolution in pixels per meter.

*biClrUsed*

Specifies the number of color indexes in the color table that are actually used by the bitmap. If this value is zero, the bitmap uses the maximum number of colors corresponding to the value of the *biBitCount* member for the compression mode specified by *biCompression*.

### *biClrImportant*

It specifies the number of color indexes that are required for displaying the bitmap. If this value is zero, all colors are required.

Phantom specific:

If *biBitCount* is 16 or 48 *biClrImportant* should be the maximum sample value + 1, which is the maximum number of sample levels.

*biClrImportant* has to be 16384, 4096 or 1024, for images having a bit depth of 14, 12 or 10 bits per sample.

### 3.6.2. IH

```
typedef struct tagIH
{
    DWORD    biSize;
    LONG     biWidth;
    LONG     biHeight;
    WORD     biPlanes;
    WORD     biBitCount;
    DWORD    biCompression;
    DWORD    biSizeImage;
    LONG     biXPelsPerMeter;
    LONG     biYPelsPerMeter;
    DWORD    biClrUsed;
    DWORD    biClrImportant;
    int      BlackLevel;
    int      WhiteLevel;
} IH, *PIH;
```

Header: PhInt.h

*biSize*

*biWidth*

*biHeight*

*biPlanes*

*biBitCount*

*biCompression*

*biSizeImage*

*biXPelsPerMeter*

*biYPelsPerMeter*

*biClrUsed*

*biClrImportant*

All of the above fields have the same meaning as their corresponding fields in the **BITMAPINFOHEADER** structure.

*BlackLevel*

Black level for the image this **IH** describes.

*WhiteLevel*

White level for the image this **IH** describes.

## 3.7. Image Processing Related Data Structures

### 3.7.1. WBGAIN

```
typedef struct tagWBGAIN
{
    float R;
    float B;
}
WBGAIN, *PWBGAIN;
```

Header: PhInt.h, PhCon.h

*R*

White balance coefficient representing gain correction for red channel.  
The gain coefficient for green channel is considered to be 1.0.

*B*

White balance coefficient representing gain correction for blue channel.

**WBGAIN** is a data structure use for storing white balance gains. It is the basic data type used when referring to an image white balance.

### 3.7.2. IMFILTER

```
typedef struct tagIMFILTER
{
    int dim;
    int shifts;
    int bias;
    int Coef[5*5];
}
IMFILTER, *PIMFILTER;
```

Header: PhInt.h

*dim*

Square convolution kernel size. It can be 3 for a 3x3 kernel or 5 for a 5x5 kernel.

*shifts*

Right shifts of *Coef* (8 shifts mean divide by 256). The shift is applied after the multiply and add operations of the convolution kernel.

*bias*

Bias to add after convolution.

*Coef*

Filter coefficients represented as integer values.  
The actual value of the coefficients is *Coef[i,j]* shifted right *shifts* times. For example if *Coef[i,j]* is 1 and *shifts* is 4 the actual float value of the coefficient is  $1/16 = 0.0625$ . *Coef* can be used both for a 3x3 or a 5x5 filter.

### 3.7.3. REDUCE16TO8PARAMS

```
typedef struct tagREDUCE16TO8PARAMS
{
    float fGain16_8;
} REDUCE16TO8PARAMS, *PREDUCE16TO8PARAMS;
```

Header: PhInt.h

*fGain16\_8*

Gain coefficient to be applied when converting to 8 bits per sample.

When converting images to 8 bits per sample, the software can take the most significant 8 bits, the least significant 8 bits or anything intermediate. To provide a fine adjustment the conversion to 8 bits will be done by multiplying the input pixel value by the factor:

$fGain16\_8 * (2^8 / 2^{bitdepth})$ .

If *fGain16\_8* is 1.0, the maximum value of the input pixel is converted to 255 (this is the maximum value for 8 bits per sample data).

### 3.7.4. TONEDESC

```
typedef struct tagToneDesc
{
    int controlPointsCounter;
    float controlPoints[32 * 2];
    char label[256];
} TONEDESC, *PTONEDESC;
```

Header: PhFile.h

*controlPointsCounter*

Number of control points describing the tone.

*controlPoints*

Tone control points array: *X*<sub>0</sub>, *y*<sub>0</sub>, ..., *X*<sub>controlPointsCounter - 1</sub>, *y*<sub>controlPointsCounter - 1</sub>,  
All values are in [0..1] range.

*label*

Optional string label describing this tone.

**TONEDESC** is the basic data type used for describing or referring to a tone.

### 3.7.5. CMDESC

```
typedef struct tagCmDesc
{
    float matrix[9];
    char label[256];
} CMDESC, *PCMDESC;
```

Header: PhFile.h

*matrix*

A 3x3 RGB color matrix.

*label*

Optional color matrix label.

**CMDESC** is the basic data type used for describing or referring to a color matrix.

## 4. User Defined Functions

### 4.1. PROGRESSCALLBACK

```
typedef BOOL (WINAPI *PROGRESSCALLBACK)(UINT CN, UINT Percent);
```

#### Parameters

*CN* [in]

Camera number.

*Percent* [in]

Percentage of work done.

#### Return Value

Value	Meaning
TRUE	Continue
FALSE	Abort

#### Remarks

PROGRESSCALLBACK depicts a callback function type a user can decide to implement and further provide as an input parameter for different functions that take a long time to complete.

The time consuming operation usually refers to a camera. That is why the camera number is one of the callback function parameters.

#### Required Header

PhCon.h

### 4.2. PROGRESSCB

```
typedef BOOL (WINAPI* PROGRESSCB)(UINT Percent, CINEHANDLE hC);
```

#### Parameters

*Percent* [in]

Percentage of work done.

*hC* [in]

Handle to the cine for which some work is in progress.

#### Return Value

Value	Meaning
TRUE	Continue saving
FALSE	Abort current save operation

#### Remarks

PROGRESSCB depicts a callback function type a user can decide to implement and further provide as an input parameter for different functions that take a long time to complete.

The time consuming operation usually refers to a cine. That is why the cine handle is one of the callback function parameters.

#### Required Header

PhFile.h

### 4.3. DEMOSAICINGFUNCPTR

```
typedef void (*DEMOSAICINGFUNCPTR) (
    PBITMAPINFOHEADER pBMIH,
    PBYTE pPixel,
    UINT CFA,
    UINT Algorithm
);
```

#### Parameters

*pBMIH* [in, out]

Input: header describing input image data.

Output: header describing output image data.

*pPixel* [in, out]

Input: Image data buffer to be color interpolated.

Output: Buffer to store color interpolated image.

*CFA* [in]

Color filter array of the camera sensor the image was recorded with.

Phantom SDK will call your function with one of the following values for the CFA parameter:

Value	Meaning
CFA_BAYER	GB RG
CFA_BAYERFLIP	RG GB
CFA_VRI	GBRG RGGB
CFA_VRIV6	BGGR GRBG

*Algorithm* [in]

Desired algorithm to be run for demosaicing the raw input image.

No Return Value

#### Remarks

This type of function is needed if an application prefers using its own demosaicing algorithms, different from the ones Phantom SDK implements. By calling `PhSetCineInfo` together with `GCI_DEMOSAICINGFUNCPTR` selector and a pointer to your own color interpolation function, the next image request will let the specified function take care of color interpolation.

Using `PhSetCineInfo/PhGetCineInfo` together with `GCI_INTALGO` selector allows an application to control which of its interpolation algorithms will be run at a given time.

Required Header

PhFile.h

## 5. Function Reference

### 5.1. Initialize/Finalize

#### 5.1.1. PhRegisterClientEx

```
HRESULT PhRegisterClientEx(
    HWND hWnd,
    char *pStgPath,
    PROGRESSCALLBACK pfnCallback,
    UINT PhConHeaderVer
);
```

##### Parameters

*hWnd* [in, optional]

Main application window used to receive notifications from the old Phantom Firewire cameras. It can be NULL.

*pStgPath* [in, optional]

Path where to store camera stg files.

It can be NULL. In this situation, a default path will be used.

*pfnCallback* [in, optional]

Pointer to a function that will be called periodically to indicate progress while searching for connected cameras.

It can be NULL.

*PhConHeaderVer* [in]

Value indicating the PhCon version the caller is using. It provides information about the version of the Phantom SDK data structures the caller application is built with.

Use the constant PHCONHEADERVERSION from PhCon.h for this parameter.

It is useful for backward compatibility.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

This function registers an application as a client of the Phantom DLLs. Call `PhRegisterClientEx` before any other function (except for `PhSetDllsOption`) whether working with Phantom cameras or with cine files.

##### Required Header

PhCon.h

#### 5.1.2. PhUnregisterClient

```
HRESULT PhUnregisterClient(HWND hWnd);
```

##### Parameters

*hWnd* [in]

Use the same value for *hWnd* as in `PhRegisterClientEx()`. It can be NULL.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

The function unregisters an application previously registered as a client of the Phantom DLLs. Call this function to ensure proper resource deallocation when the Phantom DLLs are not needed anymore.

### Required Header

`PhCon.h`

## 5.2. Create & Destroy Cine Objects

### 5.2.1. PhGetCineLive

```
HRESULT PhGetCineLive(INT CN, PCINEHANDLE phC);
```

#### Parameters

*CN* [in]

Camera number the requested live cine belongs to.

*phC* [out]

Pointer where to store live cine handle.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Returns the cinehandle associated with the live cine for the given camera. The live cine main use is reading live images from a certain camera. However, it can be used like any other cine as well; for example, for setting images processing parameters.

Please note that the live cinehandle is automatically created and automatically destroyed. `PhGetCineLive` is used just for obtaining this already allocated and initialized cinehandle.

`PhDestroyCine` must not be called for the live cine.

### Required Header

`PhFile.h`

### 5.2.2. PhNewCineFromCamera

```
HRESULT PhNewCineFromCamera(INT CN, INT CineNr, PCINEHANDLE phC);
```

#### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

RAM stored cines are indexed beginning with 1.

Flash cines are indexed beginning with value returned by the function `PhFirstFlashCine`.

*phC* [out]

Pointer where to store camera cine handle.

#### Return Value



The function returns `ERR_Ok` in case of success or `ERR_NOCINE` if there is no triggered cine in the camera. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function creates a cine handle that corresponds to a recorded camera cine. Call `PhDestroyCine` after finishing working with this cine.

#### Required Header

`PhFile.h`

#### 5.2.3. PhNewCineFromFile

```
HRESULT PhNewCineFromFile(PSTR szFN, PCINEHANDLE phC);
```

#### Parameters

*szFN* [in]

Cine file name.

*phC* [out]

Pointer where to store file cine handle.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function creates a cine handle that corresponds to a cine file. Call `PhDestroyCine` after finishing working with this cine.

#### Required Header

`PhFile.h`

#### 5.2.4. PhDuplicateCine

```
HRESULT PhDuplicateCine(PCINEHANDLE phCDest, CINEHANDLE hCSrc);
```

#### Parameters

*phCDest* [out]

Pointer or to store a new cine handle created by duplicating *hCSrc*.

*hCSrc* [in]

Source cine to be duplicated.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It creates a duplicate of the input cine.

It cannot be used for live cines.

Use this function when needing to run in parallel different operations for the same cine (e.g. reading images for display and reading images for saving to file) if cine parameters have changed from the moment the cine was created. Since some parameters are stored in internal buffers only, a `PhNewCineFromCamera` or `PhNewCineFromFile` call will also create a new cine, but it will not be aware of the parameters previously changed. This can happen especially for file cines where an image processing parameter

change is stored only in internal buffers and written to file just when the cine is saved.

Call `PhDestroyCine` after finishing working with the duplicate cine.

### Required Header

`PhFile.h`

### 5.2.5. PhGetOpenCineName

```
BOOL PhGetOpenCineName(PSTR szFileName, UINT Options);
```

#### Parameters

*szFileName* [out]

String where to store selected file name.

*Options* [in]

Available options are:

Value	Meaning
0	The user can select only one file name.
OFN_MULTISELECT	The user can select one or more file names.

#### Return Value

Value	Meaning
TRUE	The operation has succeeded. The user has ended the dialog with OK.
FALSE	The operation has been canceled. The user has ended the dialog with CANCEL.

#### Remarks

This function displays a dialog box allowing the user to select an existing cine file. The file name is returned in the parameter.

### Required Header

`PhFile.h`

### 5.2.6. PhDestroyCine

```
HRESULT PhDestroyCine(CINEHANDLE hC);
```

#### Parameters

*hC* [in]

Cine handle.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function should be called when finishing working with a specific cine. It frees any allocated resources the cine needed.

### Required Header

`PhFile.h`

## 5.3. Read Cine Images

### 5.3.1. PhGetCineImage

```
HRESULT PhGetCineImage(
    CINEHANDLE hC,
    PIMRANGE pRng,
    PBYTE pPixel,
    UINT BufferSize,
    PIH pIH
);
```

#### Parameters

*hC* [in]

Cine to read images from.

*pRng* [in, optional]

Requested image range.

For live cines, this parameter is NULL.

*pPixel* [out]

Buffer where to store requested image data.

*BufferSize* [in]

Buffer size in bytes.

If the input buffer is too small for reading and processing requested images the operation is aborted and `ERR_ImgBufferTooSmall` is returned.

*pIH* [out]

Image header describing buffer image data.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Reads and processes images from the input cine and stores them in *pPixel* buffer. The function can be used for any type of cine: live cine, camera cine or file cine.

Requested images are read from the source cine and then processed according to current use case and corresponding parameters. See `PhSetUseCase`, `PhGetUseCase` for more details.

The image header, *pIH*, fully describes the final image format.

Use *pRng* to specify what images you are interested in. For live cines, this parameter can be NULL, in which case the function will return one live image.

The pixel buffer, *pPixel*, has to be allocated by the caller. It is used for storing images all throughout the processing pipeline. For this reason, it must be large enough to accommodate all the images no matter the processing phase. To obtain the maximum memory size necessary for processing one image call `PhGetCineInfo` with the `GCI_MAXIMGSIZE` selector.

No more than one `PhGetCineImage` call per `CINEHANDLE` object must exist at any given time.

Use `PhIhtToBitmapInfoHeader` function if the image is used in a context where a `BITMAPINFOHEADER` is requested.

#### Required Header

`PhFile.h`

### 5.3.2. PhSetUseCase

```
HRESULT PhSetUseCase(CINEHANDLE hC, int CineUseCaseID);
```

#### Parameters

*hC* [in]

Cine to set the use case for.

*CineUseCaseID* [in]

Desired use case. Possible values are:

Value	Meaning
UC_VIEW	Images requested from this cine are going to be displayed/played.
UC_SAVE	Images requested from this cine are going to be written to file.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

A cine use case tells how this cine will be used and it mainly refers to cine images use. One can request cine images in order to save them to file or it can request cine images in order to display them in a video player application. Different processing steps and different optimizations apply in each of these situations. The use case represents a simple way of hiding all of these details.

By default, each cine has its use case set to `UC_VIEW`.

Please note that a cine (a `CINEHANDLE` object) cannot have more than one use case set at a time. When needing to display and save images from the same cine in parallel create one `CINEHANDLE` for each thread and set a corresponding use case for each of them. Check `PhNewCineFromCamera`, `PhNewCineFromFile` and `PhDuplicateCine` functions.

#### Required Header

`PhFile.h`

### 5.3.3. PhGetUseCase

```
HRESULT PhGetUseCase(CINEHANDLE hC, int* pCineUseCaseID);
```

#### Parameters

*hC* [in]

Cine to get the use case from.

*pCineUseCaseID* [out]

Pointer where to store the use case.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Returns the current use case for the given cine.

#### Required Header

`PhFile.h`

## 5.4. Save Cine Images to File

### 5.4.1. PhGetSaveCineName

```
BOOL PhGetSaveCineName(CINEHANDLE hC);
```

#### Parameters

*hC* [in, out]

Cine to be saved.

#### Return Value

Value	Meaning
TRUE	The operation has succeeded. The user has ended the dialog with OK.
FALSE	The operation has been canceled. The user has ended the dialog with CANCEL.

#### Remarks

This function displays a browse dialog window allowing a user to specify various cine saving options: the user can select destination file name and file path, destination format and format options or the image range to be saved. All information about the user choices is stored inside the cine. Cine saving options are also accessible by using PhSetCineInfo or PhGetCineInfo functions with corresponding selectors as described in Appendix 1.

Check also File Naming Convention chapter for special file name saving options.

PhWriteCineFile or PhWriteCineFileAsync can further be called for saving cine images to file according to previously configured options.

#### Required Header

PhFile.h

### 5.4.2. PhGetSaveCinePath

```
BOOL PhGetSaveCinePath(CINEHANDLE hC);
```

#### Parameters

*hC* [in, out]

Cine to be saved.

#### Return Value

Value	Meaning
TRUE	The operation has succeeded. The user has ended the dialog with OK.
FALSE	The operation has been canceled. The user has ended the dialog with CANCEL.

#### Remarks

This function displays a browse dialog window allowing a user to specify various options for saving a group of one or more cines: destination path, destination file naming parameters, destination format and format options. All information about the user choices is stored inside the input cine.

The caller has to further take care the same saving options are set for all the others cines in the group. This can be accomplished by calling

PhGetCineInfo and PhSetCineInfo functions with corresponding selectors as described in Appendix 1. PhWriteCineFile or PhWriteCineFileAsync can then be called for each cine in the input cine group. Check also File Naming Convention chapter for special file name saving options.

#### Required Header

PhFile.h

#### 5.4.3. PhWriteCineFile

```
HRESULT PhWriteCineFile(CINEHANDLE hC, PROGRESSCB pfnCallback);
```

##### Parameters

*hC* [in]

Cine to be saved.

*pfnCallback* [in, optional]

Pointer to the progress callback function which shows the amount of work done.

It can be NULL.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

This function saves the input cine to file. The progress callback function will provide progress updates periodically.

Before calling this function, make sure cine saving options are properly configured by calling either `PhGetSaveCineName` or `PhSetCineInfo` with corresponding selectors as described in Appendix 1.

#### Required Header

PhFile.h

#### 5.4.4. PhWriteCineFileAsync

```
HRESULT PhWriteCineFileAsync(CINEHANDLE hC, PROGRESSCB pfnCallback);
```

##### Parameters

*hC* [in]

Cine to be saved.

*pfnCallback* [in, optional]

Pointer to the progress callback function which shows the amount of work done.

It can be NULL.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

Asynchronous function used for saving a cine to file.

Asynchronous equivalent of `PhWriteCineFile`.

It starts a thread responsible for writing the file and immediately returns.

When the write operation finishes either because an error was encountered or because all data has been written to file, the callback function is called with the percent parameter set to 100. To differentiate between the two possible situations, the caller should call `PhGetCineInfo` with `GCI_WRITEERR` code selector to check if any error was encountered.

At most one save operation per cine can be active at any given time. Call `PhStopWriteCineFileAsync` to abort an asynchronous save operation.

#### Required Header

`PhFile.h`

#### 5.4.5. `PhStopWriteCineFileAsync`

```
HRESULT PhStopWriteCineFileAsync(CINEHANDLE hC);
```

##### Parameters

*hC* [in]

Cine handle.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

Stops an asynchronous save operation for the specified cine.

#### Required Header

`PhFile.h`

#### 5.4.6. `PhGetWriteCineFileProgress`

```
HRESULT PhGetWriteCineFileProgress(
    CINEHANDLE hC,
    PUINT pProgressPercent
);
```

##### Parameters

*hC* [in]

Cine handle.

*pProgressPercent* [out]

Pointer where to store current progress percent.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

Returns current save progress for the input cine.

It can be called for either an asynchronous or a synchronous write cine operation. However, please note that for a synchronous save operation the thread the save operation has been launched on is expected to be different from the thread this function is called from.

## Required Header

PhFile.h

### 5.4.7. PhGetWriteCineFileProgressAll

```
HRESULT PhGetWriteCineFileProgressAll(
    PSAVECINEINFO pSaveCineInfo,
    PUINT pCount
);
```

#### Parameters

*pSaveCineInfo* [out]

Preallocated array of `SAVECINEINFO` elements to be filled with information about all currently active write cine file operations.

*pCount* [in, out]

Input: *pSaveCineInfo* size.

Output: number of filled elements in *pSaveCineInfo* array.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Returns progress information about all synchronous and asynchronous currently active write cine file operations.

## Required Header

PhFile.h

## 5.5. Cine Parameters

### 5.5.1. PhSetCineInfo

```
HRESULT PhSetCineInfo(CINEHANDLE hC, UINT Selector, PVOID pVal);
```

#### Parameters

*hC* [in]

The cine the information is set for.

*Selector* [in]

Numerical value identifying the parameter to be set.

*pVal* [in, out]

Input: pointer containing the value to be set.

Output: pointer containing accepted value.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It is the main function used for setting cine parameters such as image processing parameters or saving options. It requires as input data a cine handle, a numerical value indicating what parameter is to be modified and the actual value to be set.



Please note that the last parameter must be passed by reference! The function will return back to the caller the value actually accepted for the specified cine parameter.  
Check Appendix 1 for a complete list of selectors `PhSetCineInfo` accepts.

#### Required Header

`PhFile.h`

#### 5.5.2. PhGetCineInfo

```
HRESULT PhGetCineInfo(CINEHANDLE hC, UINT Selector, PVOID pVal);
```

##### Parameters

*hC* [in]

The cine the information is read from.

*Selector* [in]

Numerical value identifying the parameter to be read.

*pVal* [out]

Pointer to the buffer in which the value for the requested property is to be returned.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

This function reads the cine parameter specified by *Selector*. Check Appendix 1 for a complete list of selectors `PhGetCineInfo` accepts.

#### Required Header

`PhFile.h`

#### 5.5.3. PhCineGet

```
HRESULT PhCineGet(UINT CN, INT CineNr, UINT Selector, PVOID pVal);
```

##### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

*Selector* [in]

Numerical selector identifying the parameter to be read.

*pVal* [out]

Pointer where to store requested information.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

This function gets requested information for the specified camera cine. The cine is identified based on a cine number and the camera number it belongs to.

Check Appendix 4 for possible selector values.

**Important!** This is a low level function. Use it only if neither `PhGetCineInfo` nor `PhGetCineParams` can provide needed information. There is no interaction between this function and `CINEHANDLE` objects!

#### Required Header

PhCon.h

#### 5.5.4. PhCineSet

```
HRESULT PhCineSet(UINT CN, INT CineNr, UINT Selector, PVOID pVal);
```

##### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

*Selector* [in]

Numerical selector identifying the parameter to be set.

*pVal* [in, out]

Input: value to be set.

Output: accepted value for the cine parameter.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

Sets requested information for the specified camera cine. The cine is identified based on a cine number and the camera number it belongs to. Please note that the last parameter must be passed by reference. The function will return back to the caller accepted value for the specified cine parameter.

Check Appendix 4 for possible selector values.

**Important!** This is a low level function. Use it only if neither `PhSetCineInfo` nor `PhSetCineParams`/`PhSetSingleCineParams` can provide needed functionality. There is no interaction between this function and `CINEHANDLE` objects!

#### Required Header

PhCon.h

#### 5.5.5. PhGetCineAuxData

```
HRESULT PhGetCineAuxData(
    CINEHANDLE hC,
    INT ImageNumber,
    UINT SelectionCode,
    UINT DataSize,
    PVOID pData
);
```

##### Parameters

*hC* [in]

Cine handle.

*ImageNumber* [in]

The number of the cine image.

*SelectionCode* [in]

Possible values are:

Value	Meaning
GAD_TIME	Acquisition time.
GAD_EXPOSURE	Exposure duration.
GAD_RANGE1	Range data.
GAD_BINSIG	Binary acquired signals.
GAD_ANASIG	Analog acquired signals.
GAD_SMPTE TC	SMPTE time code
GAD_SMPTE TCU	SMPTE time code unpacked

*DataSize* [in]

*pData* buffer size. It depends on the requested auxiliary data:

Requested Data	Size
GAD_TIME	sizeof(TIME64)
GAD_EXPOSURE	sizeof(UINT)
GAD_RANGE1	128 bits
GAD_BINSIG	(binary channel count + 7) / 8 x (acquired samples per image)
GAD_ANASIG	2 x (acquired samples per image) x (analog channel count)
GAD_SMPTE TC	sizeof(TC)
GAD_SMPTE TCU	sizeof(TCU)

*pData* [out]

Pointer where to store requested data.

### Return Value

The function returns `ERR_Ok` in case of success, `ERR_BadCine` if cine handle was not initialized or `ERR_AuxDataNotFound` if auxiliary data was not found. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

This function returns in *pData* requested auxiliary data for a cine.

### Required Header

PhFile.h

### 5.5.6. PhMeasureCineWB

```
HRESULT PhMeasureCineWB(
    CINEHANDLE hC,
    PPOINT pP,
    int SquareSide,
    PIMRANGE pRng,
    PWBGAIN pWB,
    PUINT pSatCnt
);
```

### Parameters

*hC* [in]

Input cine. It can be a camera or a file cine.

*pP* [in]

The point around which white balance is to be computed.

*SquareSide* [in]

Side length of the square area used to compute white balance.

*pRng* [in]

Image range to analyze when computing white balance coefficients.

*pWB* [out, optional]

Address where to store corrective gains for red and blue.

This parameter can be NULL if not needed.

*pSatCnt* [out, optional]

Address where to store saturated pixel count.

This parameter can be NULL if not needed.

### Return Value

The function returns `ERR_Ok` in case of success or `ERR_BadParam` if the input cine is not a raw cine. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

This function reads images from the input cine as specified by the *pRng* parameter and:

- if *pWB* is not NULL, it computes white balance gains
- if *pSatCnt* is not NULL, it computes average saturated pixel count

### Required Header

`PhFile.h`

### 5.5.7. PhPrintTime

```
HRESULT PhPrintTime(
    CINEHANDLE hC,
    INT ImNo,
    UINT Options,
    PSTR szTime
);
```

### Parameters

*hC* [in]

The cine for which time information is requested.

*ImNo* [in]

The cine image number the time is requested for.

*Options* [in, optional]

Numerical value specifying time format options. Possible values are:

Value	Meaning
<code>PPT_FULL</code>	Full Time
<code>PPT_DATE_ONLY</code>	Date only
<code>PPT_TIME_ONLY</code>	Time only
<code>PPT_FRAC_ONLY</code>	Fractions of second only
<code>PPT_ATTRIB_ONLY</code>	Attributes only

*szTime* [out]

Char buffer to be filled with time data according to the specified format.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

Prints the time for the given cine image according to current cine time format and specified format options. To control the time format, use `PhSetCineInfo/PhGetCineInfo` and the `GCI_TIMEFORMAT` selector.

Time format options are supported only for `TF_LT` and `TF_UT` time formats. For `TF_SMPTE` format, the time is always printed in full.

Here it is how the buffer can look like:

Time Format	Time Format Options	szTime
TF_UT	PPT_FULL	Fri May 20 2011 08:09:07.344 109.00 utc
TF_UT	PPT_DATE_ONLY	Fri May 20 2011
TF_UT	PPT_TIME_ONLY	08:09:07
TF_UT	PPT_FRAC_ONLY	.344 109.00
TF_UT	PPT_ATTRIB_ONLY	utc
TF_LT	PPT_FULL	Fri May 20 2011 11:07:23.344 109.00
TF_LT	PPT_DATE_ONLY	Fri May 20 2011
TF_LT	PPT_TIME_ONLY	11:07:23
TF_LT	PPT_FRAC_ONLY	.344 109.00
TF_LT	PPT_ATTRIB_ONLY	
TF_SMPTE	-	08:29:33:09 TC

## Required Header

PhFile.h

### 5.5.8. PhGetCineTriggerTime

```
HRESULT PhGetCineTriggerTime(CINEHANDLE hC, PTIME64 pTriggerTime);
```

#### Parameters

*hC* [in]

The cine the trigger time is requested for.

*pTriggerTime* [out]

Buffer where to store the trigger time.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Gets the moment of the trigger for the input cine as a `TIME64` value.

## Required Header

PhFile.h

## 5.6. Acquisition Parameters

### 5.6.1. PhGetCineParams

```
HRESULT PhGetCineParams(
    UINT CN,
    INT CineNr,
    PACQUIPARAMS pParams,
    PBITMAPINFO pBMI
);
```

#### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

*pParams* [out]

Pointer to an `ACQUIPARAMS` structure to be filled by this function.

*pBMI* [out]

Pointer to a `BITMAPINFO` structure to be filled by this function.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

This function is responsible for reading the acquisition parameters for a cine identified based on a cine number and a camera number for the camera the cine belongs to. It can be used both before and after the acquisition ends.

The function also fills an image header. Please allocate space for the `BITMAPINFOHEADER` and for the gray palette, `256 * sizeof(RGBQUAD)`.

Please note that this header can be different from the image header returned by `PhGetCineImage` depending on the image request parameters.

### Required Header

`PhCon.h`

### 5.6.2. PhSetCineParams

```
HRESULT PhSetCineParams(UINT CN, INT CineNr, PACQUIPARAMS pParams);
```

### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

*pParams* [in, out]

Input: buffer containing parameters to be set.

Output: buffer containing accepted acquisition parameters.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It configures the camera according to the acquisition parameters specified by *pParams* parameter. On return, this buffer contains validated and accepted values.

If the camera memory is a single partition one, use `PhSetSingleCineParams` function instead.

### Required Header

`PhCon.h`

### 5.6.3. PhSetSingleCineParams

```
HRESULT PhSetSingleCineParams(UINT CN, PACQUIPARAMS pParams);
```

### Parameters

*CN* [in]

Camera number.

*pParams* [in, output]

Input: buffer containing parameters to be set.

Output: buffer containing accepted acquisition parameters.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It configures the camera according to the acquisition parameters specified by *pParams* parameter. On return, this buffer contains validated and accepted values.

This function is to be used for single partition mode only.

### Required Header

PhCon.h

#### 5.6.4. PhGetBitDepths

```
HRESULT PhGetBitDepths(UINT CN, PUINT pCnt, PUINT pBitDepths);
```

### Parameters

*CN* [in]

Camera number.

*pCnt* [out]

Count of values in *pBitDepths*.

*pBitDepths* [out]

Array to be filled with possible recording bit depths for this camera.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It provides a list of the possible pixel bit depths when recording in the camera FBM. Example: the Phantom v9.1 camera can record at 8, 10, 12, 14, 16 bits per pixel. Larger bit depths may provide better image or better sensitivity but reduce the number of images that can be recorded in the camera memory.

The bit depths are stored in the array in increasing order.

### Required Header

PhCon.h

#### 5.6.5. PhGetResolutions

```
HRESULT PhGetResolutions(
    UINT CN,
    PPOINT pRes,
    PUINT pCnt,
    PBOOL pCAR,
```

```
    PUINT pADCBits
);
```

## Parameters

*CN* [in]

Camera number.

*pRes* [out, optional]

Array of `POINT` elements to be filled with resolutions available for this camera.

It can be `NULL` if this information is not needed, but in this case *pCnt* must be `NULL` as well.

*pCnt* [in, out, optional]

Input: Pointer containing *pRes* number of allocated elements.

Output: Pointer containing *pRes* number of filled elements.

It can be `NULL` if this information is not needed, but in this case *pRes* must be `NULL` as well.

*pCAR* [in, out]

At input, if *pCAR* is `NULL` the result in *pRes* will be a list of usual resolutions even if the camera supports Continuous Adjustable Resolution (CAR).

At output, if *pCAR* is not `NULL`:

\**pCAR* will be `TRUE` if the camera supports CAR.

In this case the results in *pRes* are:

Value	Meaning
<i>pRes</i> [0]	Maximum Resolution (compatible with the old style).
<i>pRes</i> [1]	Minimum Resolution.
<i>pRes</i> [2]	Increments.

Available (X, Y) resolutions can further be generated following these guidelines:

```
X = MinX + i * IncX;    i = 0..M
Y = MinY + j * IncY;    j = 0..N
```

where:

```
MaxX = pRes[0].x
MaxY = pRes[0].y
MinX = pRes[1].x
MinY = pRes[1].y
IncX = pRes[2].x
IncY = pRes[2].y
M = (MaxX - MinX) / IncX
N = (MaxY - MinY) / IncY
```

*pADCBits* [out, optional]

Pointer where to store the maximum bit depth available for this camera.

This is actually the analog to digital converters resolution.

It can be `NULL` if this information is not needed.

## Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

## Remarks

It gets a list of image dimensions available for this camera, as well as the maximum pixel depth.

If one resolution is requested and several are available, the maximum camera resolution is returned in *pRes*.



## Required Header

PhCon.h

### 5.6.6. PhGetFrameRateRange

```
HRESULT PhGetFrameRateRange(
    UINT CN,
    POINT Res,
    UINT *pMinFrameRate,
    UINT *pMaxFrameRate
);
```

#### Parameters

*CN* [in]

Camera number.

*Res* [in]

Image reference dimension.

*pMinFrameRate* [out]

Pointer where to store the minimum frame rate.

*pMaxFrameRate* [out]

Pointer where to store the maximum frame rate.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Given a reference resolution for a camera, this routine finds the minimum and the maximum available frame rates.

## Required Header

PhCon.h

### 5.6.7. PhGetExactFrameRate

```
HRESULT PhGetExactFrameRate(
    UINT CamVer,
    UINT SyncMode,
    UINT RequestedFrameRate,
    double *pExactFrameRate
);
```

#### Parameters

*CamVer* [in]

Camera version.

*SyncMode* [in]

Source of image acquisition pulses.

It has the same value as the *SyncImaging* field of an `ACQUIPARAMS` structure containing current acquisition parameters.

*RequestedFrameRate* [in]

Requested frame rate.

*pExactFrameRate* [out]

The exact frame rate the camera can achieve, close to the requested one.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

Given a requested frame rate, this function computes the exact frame rate the camera is able to accomplish. The period corresponding to this frame rate has to be an exact multiple of the internal clock period of the camera. This routine does not check if the resulted rate is in the interval (minimum frame rate, maximum frame rate). Call `PhGetFrameRateRange` to validate returned frame rate.

### Required Header

PhCon.h

### 5.6.8. PhGetExposureRange

```
HRESULT PhGetExposureRange(
    UINT CN,
    UINT FrameRate,
    UINT *pMinExposure,
    UINT *pMaxExposure
);
```

### Parameters

*CN* [in]

Camera number.

*FrameRate* [in]

Frame rate.

*pMinExposure* [out]

Pointer where to store minimum exposure (nanoseconds).

*pMaxExposure* [out]

Pointer where to store maximum exposure(nanoseconds).

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

This function computes the limits of the exposure parameter for the specified frame rate.

### Required Header

PhCon.h

### 5.6.9. PhParamsChanged

```
HRESULT PhParamsChanged(UINT CN, PBOOL pChanged);
```

### Parameters

*CN* [in]

Camera number.

*pChanged* [in]

If `TRUE`, changes have been made to the camera's acquisition parameters.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Use this function to determine whether another computer, remote control, or on-camera control unit changed the acquisition or image processing parameters of the camera. If *pChanged* is `TRUE`, changes have been made and the new settings must be read from camera.

#### Required Header

PhCon.h

## 5.7. Camera Pool

### 5.7.1. PhGetCameraCount

```
HRESULT PhGetCameraCount (UINT *pCnt);
```

#### Parameters

*pCnt* [out]

Pointer where to store camera count.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It gets the count of Phantom cameras connected to the computer. If there are no connected cameras (or if you access a camera number above or equal to this value) the library switches to a simulated mode that allows you to test or to use the software without a camera connected.

#### Required Header

PhCon.h

### 5.7.2. PhConfigPoolUpdate

```
HRESULT PhCheckCameraPool (UINT Period);
```

#### Parameters

*Period* [in]

Interval between camera enumerations in milliseconds.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

`PhConfigPoolUpdate()` will enable automatic search for new or revived cameras. It will create a separate thread for enumeration of the cameras in Phantom SDK. *Period* will define the interval between two broadcasts. If *Period* is 0, the function will stop the automatic enumeration (if enabled) and will do a single re-enumeration.

The periodic or single enumerations mentioned above will do both identifying there is a change in the camera pool and will activate that change (revive an Offline camera or add a new camera to the end of the pool).

## Required Header

PhCon.h

### 5.7.3. PhGetAllIpCameras

```
HRESULT PhGetAllIpCameras(PUINT pCnt, PCAMERAID pCam);
```

#### Parameters

*pCnt* [out]

Address where to store discovered camera count.

*pCam* [out]

Buffer to be filled with CAMERAID elements describing discovered cameras.

#### Return Value

The function returns ERR\_Ok in case of success or an error code, otherwise. To get the text of the error message call PhGetErrorMessage.

#### Remarks

It provides an option for getting information about connected cameras and can be called in parallel with other camera operations. It does not affect or interact with current camera pool.

## Required Header

PhCon.h

### 5.7.4. PhOffline

```
BOOL PhOffline(UINT CN);
```

#### Parameters

*CN* [in]

Camera number

#### Return Value

The function returns TRUE when the connection to the camera was lost and FALSE if the Ethernet camera work normally.

#### Remarks

When an attempted access to a camera returns a timeout error, the camera is not removed from the pool but its status is changed to "Offline". If the connection to the camera is working again and PhConfigPoolUpdate() activated the enumeration thread, the camera is revived and can be accessed.

## Required Header

PhCon.h

### 5.7.5. PhMakeAllIpVisible

```
HRESULT PhMakeAllIpVisible();
```

#### No Parameters

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Empties both the *Visible* and the *Ignored* camera lists. This way, any camera will become visible.

#### Required Header

PhCon.h

#### 5.7.6. PhGetVisibleIp

```
HRESULT PhGetVisibleIp(PUINT pCnt, PUINT pIpAddress);
```

#### Parameters

*pCnt* [out]

Count of visible camera IP addresses.  
Limited to `MAXCAMERACNT` (see PhCon.h).

*pIpAddress* [out]

Array of visible camera IP addresses.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Use this function to obtain an array with IP addresses for all the cameras in the visible or "Connect Only To" list. The visible list contains the cameras you wish to be connected to, no matter if they are or are not on the network of your computer at that moment.

If the visible list contains at least one camera, the computer connects only to cameras in the visible list and does not connect to cameras in the ignored list.

If you include a camera in both the visible and ignored lists, the visible list will have priority so the camera will be connected to your computer.

#### Required Header

PhCon.h

#### 5.7.7. PhAddVisibleIp

```
HRESULT PhAddVisibleIp(UINT IpAddress);
```

#### Parameters

*IpAddress* [in]

Camera IP address.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It adds a camera IP address to the visible list. The visible list cannot have more than `MAXCAMERACNT` (see PhCon.h) elements.

#### Required Header

PhCon.h

#### 5.7.8. PhRemoveVisibleIp

```
HRESULT PhRemoveVisibleIp(UINT IpAddress);
```

##### Parameters

*IpAddress* [in]  
Camera IP address.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

It removes a camera IP address from the visible list.

##### Required Header

PhCon.h

#### 5.7.9. PhDisableVisible

```
HRESULT PhDisableVisible(BOOL bDisabled);
```

##### Parameters

*bDisabled* [in]  
TRUE/FALSE enable or disable visible list.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

Enable or disable the visible list. It allows the list to remain persistent when you want to switch temporary to "See all cameras" mode.

##### Required Header

PhCon.h

#### 5.7.10. PhIsVisibleDisabled

```
HRESULT PhIsVisibleDisabled(void);
```

##### Parameters

##### Return Value

The function returns the current status of the Visible list.

##### Remarks

The visible list can be active or disabled.

##### Required Header

PhCon.h

#### 5.7.11. PhFillIDInfo

```
HRESULT PhFillIDInfo(PCAMERAID pCam);
```

##### Parameters

*pCam* [in, out]

Camera identification information including IP address, name, serial, model .... IP is "in" parameter for this function, all other are "out".

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

You have to fill the IP address of the camera `pCam->IP` ( it has to be net order i.e. big endian here) and the function will fill the Serial, Name, Model, CFA, CamVer. This information is provided even if the camera is absent now but was seen in the past at least once. It can be used to populate your list of cameras with more friendly details about the camera.

##### Required Header

PhCon.h

#### 5.7.12. PhGetIgnoredIp

```
HRESULT PhGetIgnoredIp(PUINT pCnt, PUINT pIpAddress);
```

##### Parameters

*pCnt* [out]

Count of ignored camera IP addresses.  
Limited to `MAXCAMERACNT` (see PhCon.h).

*pIpAddress* [out]

Array of ignored camera IP addresses.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

Use this function to obtain an array with IP addresses for all cameras in the ignored, or "Do Not Connect To" list. The ignored list contains the cameras which you do not wish to be connected to, no matter if they are or are not on the network of your computer at that moment.

##### Required Header

PhCon.h

#### 5.7.13. PhAddIgnoredIp

```
HRESULT PhAddIgnoredIp(UINT IpAddress);
```

##### Parameters

*IpAddress* [in]

Camera IP address.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It adds a camera IP address to the ignored list. The ignored list cannot have more than `MAXCAMERACNT` (see `PhCon.h`) elements.

### Required Header

`PhCon.h`

### 5.7.14. PhRemoveIgnoredIp

```
HRESULT PhRemoveIgnoredIp(UINT IpAddress);
```

### Parameters

*IpAddress* [in]

Camera IP address.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It removes a camera IP address from the ignored list.

### Required Header

`PhCon.h`

### 5.7.15. PhAddSimulatedCamera

```
HRESULT PhAddSimulatedCamera(UINT CamVer, UINT Serial);
```

### Parameters

*CamVer* [in]

Simulated camera version.

*Serial* [in]

Desired serial for the simulated camera.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It adds a new simulated camera to the camera pool. `PhGetCameraCount` will count these cameras too.

You may later load a stg file for the camera having the specified serial and simulate the behavior and the calibrations of that camera.

### Required Header

`PhCon.h`



## 5.8. Camera Parameters

### 5.8.1. PhGet

```
HRESULT PhGet(UINT CN, UINT Selector, PVOID pVal);
```

#### Parameters

*CN* [in]

Camera number.

*Selector* [in]

Selector value indicating the information to be retrieved.

*pVal* [out]

Pointer to the buffer in which the value for the requested property is to be returned.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function gets information about the given camera based on a selection code. For a complete list of currently available selectors, go to Appendix 3.

#### Required Header

PhCon.h

### 5.8.2. PhSet

```
HRESULT PhSet(UINT CN, UINT Selector, PVOID pVal);
```

#### Parameters

*CN* [in]

Camera number.

*Selector* [in]

Selector value indicating information to be set.

*pVal* [out]

Pointer to the buffer containing data to be set.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function performs a set operation for the given camera. For a complete list of currently available selectors, go to Appendix 3. Please note that the last parameter must be passed by reference. The function will return back to the called the value actually accepted for the specified parameter.

#### Required Header

PhCon.h

### 5.8.3. PhGet1

```
HRESULT PhGet1(UINT CN, UINT Selector, UINT Selector1, PVOID pVal);
```

### Parameters

*CN* [in]

Camera number.

*Selector* [in]

Selector value indicating the information to be retrieved.

*Selector1* [in]

Programmable port number or other information, depending on *Selector*  
Range of possible values is [1..N] where <N> is the value returned from  
PhGet using *gsPortCount* selector

*pVal* [out]

Pointer to the buffer in which the value for the requested property is to be returned.

### Return Value

The function returns *ERR\_Ok* in case of success or an error code,  
otherwise. To get the text of the error message call *PhGetErrorMessage*.

### Remarks

This function gets information about the given camera based on a  
selection code. For a complete list of currently available selectors, go to  
Appendix 5.

### Required Header

PhCon.h

#### 5.8.4. PhSet1

```
HRESULT PhSet1(UINT CN, UINT Selector, UINT Selector1, PVOID pVal);
```

### Parameters

*CN* [in]

Camera number.

*Selector* [in]

Selector value indicating information to be set.

*Selector1* [in]

Programmable port number or other information, depending on *Selector*  
Range of possible values is [1..N] where <N> is the value returned from  
PhGet using *gsPortCount* selector

*pVal* [out]

Pointer to the buffer containing data to be set.

### Return Value

The function returns *ERR\_Ok* in case of success or an error code,  
otherwise. To get the text of the error message call *PhGetErrorMessage*.

### Remarks

This function performs a set operation for the given camera.  
For a complete list of currently available selectors, go to Appendix 5.  
Please note that the last parameter must be passed by reference. The  
function will return back to the called the value actually accepted for the  
specified parameter.

### Required Header

PhCon.h

### 5.8.5. PhGet2

```
HRESULT PhGet2(
    UINT CN,
    UINT Selector,
    UINT Selector1,
    UINT Selector2
    PVOID pVal
);
```

#### Parameters

*CN* [in]

Camera number.

*Selector* [in]

Selector value indicating the information to be retrieved.

Name	Meaning	Type	Get	Set
gsSigName	Signal name, given the port number and signal number	Char[ <b>MAXSTDSTRSZ</b> ]	✓	✗

*Selector1* [in]

Programmable port number

Range of possible values is [1..N] where <N> is the value returned from PhGet using gsPortCount selector

*Selector2* [in]

Programmable port signal index

Range of possible values is [0..M-1] where <M> is the value returned from PhGet using gsSigCount selector

*pVal* [out]

Pointer to the buffer in which the value for the requested property is to be returned.

#### Return Value

The function returns **ERR\_Ok** in case of success or an error code, otherwise. To get the text of the error message call **PhGetErrorMessage**.

#### Remarks

This function gets information about the given camera based on a selection code.

#### Required Header

PhCon.h

### 5.8.6. PhGetCameraOptions

```
HRESULT PhGetCameraOptions(UINT CN, PCAMERAOPTIONS pOptions);
```

#### Parameters

*CN* [in]

Camera number.

*pCAMERAOPTIONS* [out]

Buffer to be filled with current camera options.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function reads general camera configuration parameters and stores them in the `CAMERAOPTIONS` structure the caller provided.

#### Required Header

PhCon.h

#### 5.8.7. PhSetCameraOptions

```
HRESULT PhSetCameraOptions(UINT CN, PCAMERAOPTIONS pOptions);
```

#### Parameters

*CN* [in]

Camera number.

*pCAMERAOPTIONS* [in]

Pointer to the camera options.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It configures the camera according to parameters specified by the `CAMERAOPTIONS` structure provided.

#### Required Header

PhCon.h

#### 5.8.8. PhGetVersion

```
HRESULT PhGetVersion(UINT CN, UINT VerSel, PUINT pVersion);
```

#### Parameters

*CN* [in]

Camera number.

*VerSel* [in]

Code to select the desired information.

*pVersion* [out]

Buffer where to store requested information.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It reads version information for the input camera as indicated by the *VerSel* parameter. *VerSel* parameter can have one of the following values:

Value	Meaning
GV_CAMERA	Camera hardware version.
GV_FIRMWARE	Camera firmware version.
GV_FPGA	Camera FPGA version.
GV_PHCON	PhCon.dll version.
GV_CFA	Sensor CFA.

GV_KERNEL	Camera kernel version.
GV_MAGAZINE	Camera magazine version.
GV_FIRMWAREPACK	Camera firmware packed version. (Ph16 Cameras)
GV_HEAD_FIRMWARE	Miro N5 head firmware version.
GV_HEAD_FPGA	Miro N5 FPGA.
GV_HEAD_FIRMWAREPACK	Miro N5 firmware packed version.
GV_VSHUTTER_FIRMWARE	V Shutter firmware version.

## Required Header

PhCon.h

### 5.8.9. PhGetCameraID

```
HRESULT PhGetCameraID(UINT CN, UINT *pSerial, char *pCameraName);
```

#### Parameters

*CN* [in]

Camera number.

*pSerial* [out, optional]

Pointer where to store the serial number of the camera.

It can be NULL.

*pCameraName* [out, optional]

Buffer where to store camera name.

Maximum 15 characters.

It can be NULL.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It reads camera name and camera factory serial number.

The serial number is unique for each camera while the camera name can be changed by the user.

## Required Header

PhCon.h

### 5.8.10. PhSetCameraName

```
HRESULT PhSetCameraName(UINT CN, char *pCameraName);
```

#### Parameters

*CN*[in]

Camera number.

*pCameraName* [in]

Buffer containing camera name on maximum 15 characters.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It allows the user to change the camera's name.

## Required Header

PhCon.h

#### 5.8.11. PhFindCameraNumber

```
HRESULT PhFindCameraNumber(UINT Serial, UINT *pCN);
```

##### Parameters

*Serial* [in]

Factory serial number.

*pCN* [out]

Pointer where to store camera number.

##### Return Value

The function returns `ERR_Ok` in case of success or `ERR_SerialNotFound` if it does not find the camera with the specified serial.

To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

It scans all cameras inside the camera pool to find the one with the specified serial number.

##### Required Header

PhCon.h

#### 5.8.12. PhGetCameraModel

```
HRESULT PhGetCameraModel(UINT CamVer, PSTR pModel);
```

##### Parameters

*CamVer* [in]

Camera version.

*pModel* [out]

Buffer where to store a text description of the camera model.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

Use this function to obtain a text description corresponding to an integer camera model value. For example, when *CamVer* is 120 the result in *pModel* will be "Phantom v12".

##### Required Header

PhCon.h

#### 5.8.13. PhGetCameraTime

```
HRESULT PhGetCameraTime(UINT CN, PTIME64 pTime64);
```

##### Parameters

*CN* [in]

Camera number.

*pTime64* [out]

Buffer where to put the current time in the `TIME64` format.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It reads the camera clock.

The fractions of second are important for the information stored in bit 0 and bit 1: the presence of IRIG synchronization signal and the Event binary input.

Use the standard `ctime` function to convert it to a string of characters.

### Required Header

`PhCon.h`

#### 5.8.14. PhSetCameraTime

```
HRESULT PhSetCameraTime(UINT CN, UINT32 Time);
```

### Parameters

*CN* [in]

Camera number.

*Time* [in]

Seconds from January 1, 1970 0:0:0.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It sets the camera clock (for example from the computer clock).

If IRIG signal is present, this operation will fail.

### Required Header

`PhCon.h`

#### 5.8.15. PhMaxCineCnt

```
HRESULT PhMaxCineCnt(UINT CN);
```

### Parameters

*CN* [in]

Camera number.

### Return Value

The function returns the maximum number of cines that can be in the camera.

### Remarks

Phantom cameras have a preview cine numbered as 0 and a number of recordable cines numbered from 1 upward. Depending on camera the maximum count of recordable cines can be 63, 511. This function replace the constant `MAXCINECNT` from previous versions of the SDK.

### Required Header

`PhCon.h`

#### 5.8.16. PhFirstFlashCine

```
HRESULT PhFirstFlashCine(UINT CN);
```

##### Parameters

*CN* [in]

Camera number.

##### Return Value

The function returns the offset for the numbers of the cines available in the camera flash (magazine, SSD).

##### Remarks

The function replaces the previous constant `FIRST_FLASH_CINE` and has to be used to compute the cine numbers to access the recordings from flash.

##### Required Header

PhCon.h

#### 5.8.17. PhGetPartitions

```
HRESULT PhGetPartitions(UINT CN, PUINT pCount, PUINT pPartitionSize);
```

##### Parameters

*CN* [in]

Camera number.

*pCount* [in, out]

Input: Pointer containing *pPartitionSize* number of allocated elements.

Output: Pointer containing *pPartitionSize* number of elements filled by this function

*pPartitionSize* [out]

Array to be filled with memory size for each camera partition.

##### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

##### Remarks

This function fills the input array with memory size in MB for each camera partition.

##### Required Header

PhCon.h

#### 5.8.18. PhSetPartitions

```
HRESULT PhSetPartitions(UINT CN, UINT Count, PUINT pWeights);
```

##### Parameters

*CN* [in]

Camera number.

*Count* [in]

Partition count.

*pWeights* [in]



Array of weights of the memory size to be allocated to each cine. To get the percentage of each partition, it multiplies its weight by 100 and divides the result by the sum of all weights.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function is used for partitioning a camera memory. The caller must specify desired partition count as well as a weight for each partition size. All recorded cines are deleted when this function is called.

#### Required Header

`PhCon.h`

### 5.8.19. PhGetCineStatus

```
HRESULT PhGetCineStatus(UINT CN, PCINESTATUS pStatus);
```

#### Parameters

*CN* [in]

Camera number.

*pStatus* [out]

Array of `CINESTATUS` elements to be filled with information about each cine status. The buffer must be allocated for the maximum cine count *ij* the camera obtained using `PhMaxCineCnt()`.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It checks the status for all cines in this camera (finished, active, or triggered).

#### Required Header

`PhCon.h`

### 5.8.20. PhGetCineCount

```
HRESULT PhGetCineCount(UINT CN, PUINT pRAMCount, PUINT pNVMCount);
```

#### Parameters

*CN* [in]

Camera number.

*pRAMCount* [out]

Pointer where to store recorded cine count.

*pNVMCount* [out]

Pointer where to store NVM stored cine count.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It counts the RAM and NVM cines available for playback from the specified camera.

#### Required Header

PhCon.h

#### 5.8.21. PhMemorySize

```
HRESULT PhMemorySize(UINT CN, PUINT pDRAMSize, PUINT pNVMSize);
```

#### Parameters

*CN* [in]

Camera number.

*pDRAMSize* [out]

Pointer where to store the DRAM memory size in MB.

*pNVMSize* [out]

Pointer where to store the NonVolatile memory size in MB.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It provides the size of the camera DRAM and NonVolatile Memory.

#### Required Header

PhCon.h

### 5.9. Camera Control

#### 5.9.1. PhRecordCine

```
HRESULT PhRecordCine(UINT CN);
```

#### Parameters

*CN* [in]

Camera number.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It starts a new recording in the first partition of the camera.  
Any previously recorded cine will be overwritten. Please note that this operation may result in unwanted data loss if the initial cine is not saved.

#### Required Header

PhCon.h

#### 5.9.2. PhRecordSpecificCine

```
HRESULT PhRecordSpecificCine(UINT CN, int CineNr);
```

#### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number for which to start the recording.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It starts a new recording in the specified partition of a multicine camera. It deletes the existing cine if there is one recorded.

#### Required Header

`PhCon.h`

### 5.9.3. PhSendSoftwareTrigger

```
HRESULT PhSendSoftwareTrigger(UINT CN);
```

#### Parameters

*CN* [in]

Camera number.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It sends a software trigger to the camera. This is a handy way to finish a recording.

#### Required Header

`PhCon.h`

### 5.9.4. PhDeleteCine

```
HRESULT PhDeleteCine(UINT CN, INT CineNr);
```

#### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It deletes the specified camera cine. A new recording may use the parameters and the space of that cine.

#### Required Header

`PhCon.h`

### 5.9.5. PhVideoPlay

```
HRESULT PhVideoPlay(
    UINT CN,
    INT CineNr,
    PIMRANGE pRng,
    INT PlaySpeed
);
```

#### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number for a recorded camera cine.

*pRng* [in, optional]

Range of cine images to be played.

If this is `NULL`, the function changes only the playback speed to the value specified in the *PlaySpeed* parameter.

*PlaySpeed* [in]

Playback speed in milihertz.

Negative values means play backward:

Value	Meaning
1000	1 frame per second
< 0	Play backward.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It plays back a clip from the specified cine to the camera video out.

Parameter Values	Meaning
<code>CineNr &lt; 1</code>	The video will show preview images.
<code>pRng-&gt;Cnt == 0</code> or <code>PlaySpeed == 0</code>	Continuously play <code>pRng-&gt;First</code> image (Pause Mode). Using this feature, one can implement Step Forward / Step Backward.
<code>pRng-&gt;Cnt != 0</code>	Playback the specified range of images to the video output forward (if <code>PlaySpeed</code> is positive) or backward (if <code>PlaySpeed</code> is negative).

#### Required Header

PhCon.h

### 5.9.6. PhGetVideoFrNr

```
HRESULT PhGetVideoFrNr(UINT CN, PINT pCrtIm);
```

#### Parameters

*CN* [in]

Camera number.

*pCrtIm* [out]

Pointer where to store current image number being played back to the video output.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It gets the image number that is played back to the video out at the current moment. It can be used to update a browse slider or a numeric display in a video player.

#### Required Header

PhCon.h

#### 5.9.7. PhGetCameraErrMsg

```
HRESULT PhGetCameraErrMsg(UINT CN, char *pErrMsg);
```

#### Parameters

*CN* [in]

Camera number.

*pErrMsg* [out]

Error message.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It provides the last error message from the camera firmware. After calling the function, the error buffer of this camera is emptied.

#### Required Header

PhCon.h

### 5.10. Camera Calibration

#### 5.10.1. PhBlackReference

```
HRESULT PhBlackReference(UINT CN, PROGRESSCALLBACK pfnCallback);
```

#### Parameters

*CN* [in]

Camera number.

*pfnCallback* [in]

Progress indicator function.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function is responsible for controlling a new black reference operation. Images are acquired in the darkness and averaged in order to determine pixel offset values. These offsets are stored in internal tables, in the stg file and, depending on the camera version, inside the camera as well.

During a black reference operation some acquisition parameters may be changed, but they will be restored to their initial values after the calibration finishes.

Because images must be acquired in the dark, camera lenses must be covered while running a black calibration. This is not a problem for cameras having a mechanical shutter: PhBlackReference will take care of controlling the shutter as needed. Call PhGet function with the gsHasMechanicalShutter selector to find out if a camera has a mechanical shutter or not.

After the black reference operation finished, call PhWriteStgFlash for the new offset values to be sent to camera. This way the same calibration values will be used for video output.

Better results can be obtained by calling PhBlackReferenceCI which does the calibration without changing the resolution or other acquisition parameters.

### Required Header

PhCon.h

### 5.10.2. PhBlackReferenceCI

```
HRESULT PhBlackReferenceCI(UINT CN, PROGRESSCALLBACK pfnCallback);
```

#### Parameters

CN [in]

Camera number.

pfnCallback [in]

Progress indicator function.

#### Return Value

The function returns ERR\_Ok in case of success or an error code, otherwise. To get the text of the error message call PhGetErrorMessage.

#### Remarks

This function is responsible for controlling a new black reference operation at the current acquisition parameters, also names *Current Session Reference* or *CSR*.

Images are acquired in the darkness and averaged in order to determine pixel offset values. These offsets are to be used for the current acquisition parameters only. Any acquisition parameter change may require repeating the black reference operation.

Because images must be acquired in the dark, camera lenses must be covered while running a black calibration. This is not a problem for cameras having a mechanical shutter: PhBlackReferenceCI will take care of controlling the shutter as needed. Call PhGet function with the gsHasMechanicalShutter selector to find out if a camera has a mechanical shutter or not.

Newer cameras are able to perform themselves this operation. Offsets are stored inside the camera and automatically applied to image pixels. This way, the new calibration will immediately take effect on the video output. However, for the new offsets to take effect on older cameras' video output as well, they must be explicitly written to camera by calling PhWriteStgFlash function.

To distinguish between the two types of cameras call PhGet function with the gsSupportsOffGainCorrections selector.

## Required Header

PhCon.h

### 5.10.3. PhWriteStgFlash

```
HRESULT PhWriteStgFlash(UINT CN, PROGRESSCALLBACK pfnCallback);
```

#### Parameters

*CN* [in]

Camera number.

*pfnCallback* [in]

Pointer to the progress indicator function. This operation may take longer time for high resolution cameras.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Use this function to send current calibration information to the camera flash in order to be used for the video output.

This operation is necessary:

- after a CSR using `PhBlackReferenceCI` if the current camera does not implement offset gain corrections internally (It returns `FALSE` on calling `PhGet` function with the `gsSupportsOffGainCorrections` selector)
- after performing a black calibration using `PhBlackReference`

## Required Header

PhCon.h

### 5.11. NVM

#### 5.11.1. PhNVMGetStatus

```
HRESULT PhNVMGetStatus(
    UINT CN,
    PUINT pCineCnt,
    PUINT32 pTime,
    PUINT pFreeSp
);
```

#### Parameters

*CN* [in]

Camera number.

*pCineCnt* [in, out]

Input: *pTime* number of elements.

Output: Count of cines stored in NV memory.

*pTime* [out, optional]

Array with the trigger time in seconds for each cine.

It can be NULL.

*pFreeSp* [out]

Free space in NVM in MB.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function reads the number of NVM stored cines and, optionally, NVM free space. Also optionally, *pTime* array can be filled with each cine trigger time represented in seconds.

#### Required Header

PhCon.h

#### 5.11.2. PhNVMerase

```
HRESULT PhNVMerase(UINT CN, PROGRESSCALLBACK pfnCallback);
```

#### Parameters

*CN* [in]

Camera number.

*pfnCallback* [in]

Callback function to see the progress.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It erases all the cines from camera NVM.

#### Required Header

PhCon.h

#### 5.11.3. PhNVMContRec

```
HRESULT PhNVMContRec(
    UINT CN,
    PBOOL pSave,
    PUINT pOldValue,
    PUINT pEnable
);
```

#### Parameters

*CN* [in]

Camera number.

*pSave* [out]

If `TRUE`, this camera is saving the cine to the NVM or is automatically playing back to video so the images cannot be downloaded – that happens at Firewire cameras.

*pOldValue* [out, optional]

Buffer where to store previous status information. It is a combination of the following flags (see PhCon.h):

Value	Meaning
<code>NVCR_APV</code>	Automatic playback to video.
<code>NVCR_CONT_REC</code>	Continuous recording to NVM mode.
<code>NVCR_REC_ONCE</code>	Recording once mode.

*pEnable* [in, optional]



Pointer where a numerical value is stored. This numerical value controls the continuous recording mode. It is a combination between values described for *pOldValue*.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It sets the continuous recording mode for the nonvolatile memory. When a cine is available it is saved in the NV memory and the recording of a new cine is automatically started. This function tests if a cine save to NVM operation is in progress in a Firewire camera using the `pSave`.

An easier way of controlling these automatic actions is to configure *AutoSaveNVM*, *AutoSaveFile*, *AutoPlay* and *AutoCapture* fields of the `CAMERAOPTIONS` structure and further call `PhSetCameraOptions` function.

#### Required Header

`PhCon.h`

#### 5.11.4. PhNVMGetSaveRange

```
HRESULT PhNVMGetSaveRange(UINT CN, PIMRANGE pRng);
```

#### Parameters

*CN* [in]

Camera number.

*pRng* [out]

Image range.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It gets the current image range configured for a NVM Save or NVM continuous recording.

#### Required Header

`PhCon.h`

#### 5.11.5. PhNVMSaveRange

```
HRESULT PhNVMSaveRange(UINT CN, PIMRANGE pRng);
```

#### Parameters

*CN* [in]

Camera number.

*pRng* [in]

Image range.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

It sets the image range to be used for a NVM save or NVM continuous recording. It does not affect the range used by PhNVMSaveClip.

### Required Header

PhCon.h

### 5.11.6. PhNVMSave

```
HRESULT PhNVMSave(UINT CN, INT CineNr, PROGRESSCALLBACK pfnCallback);
```

#### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

*pfnCallback* [in]

Callback function to see the progress.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Saves the specified DRAM recorded cine to camera non-volatile memory.

### Required Header

PhCon.h

### 5.11.7. PhNVMSaveClip

```
HRESULT PhNVMSaveClip(
    UINT CN,
    int CineNr,
    PIMRANGE pRng,
    UINT Options,
    PROGRESSCALLBACK pfnCallback
);
```

#### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

*pRng* [in]

Image range to be saved.

`NULL` pointer means save full cine.

*Options* [in]

Use `GI_BPP16` to save the cine on 16 bits. If absent, the bit depth will be 8 bits.

*pfnCallback* [in]

Callback function to see the progress.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

Saves the specified DRAM recorded cine to camera non-volatile memory. This saving does not depend on the range set by `PhNVMSetSaveRange`.

### Required Header

PhCon.h

### 5.11.8. PhNVMRestore

```
HRESULT PhNVMRestore(
    UINT CN,
    UINT CineNr,
    PROGRESSCALLBACK pfnCallback
);
```

### Parameters

*CN* [in]

Camera number.

*CineNr* [in]

Cine number.

*pfnCallback* [in]

Callback function to see the progress.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It transfers the specified cine from the NVMemory to DRAM. Only for the Firewire cameras. In all Ethernet cameras NVM cines can be read the same way as FBM cines. However, for compatibility to the previous versions a „virtual restore“ is available from PhCon.

### Required Header

PhCon.h

## 5.12. DLLs Options

### 5.12.1. PhSetDllsOption

```
HRESULT PhSetDllsOption(UINT optionSelector, PVOID pValue);
```

### Parameters

*optionSelector* [in]

Value indicating the option to be set. Possible values are:

Value	Meaning
DO_IGNORECAMERAS	Ignore all cameras when a new application registers as a DLL client.

*pValue* [out]

Buffer containing the value to be set. This value can have one of the following types:

Selector Value	Type
DO_IGNORECAMERAS	INT

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Sets various user options regarding Phantom dlls use. This function must be called before `PhRegisterClientEx`.

#### Required Header

`PhCon.h`

#### 5.12.2. PhGetErrorMessage

```
HRESULT PhGetErrorMessage(int ErrNo, char *pErrText);
```

#### Parameters

*ErrNo* [in]

Error code returned by Phantom functions.

*pErrText* [out]

A text that describes the error.

Allocate this buffer for `MAXERRMESS` (see `PhCon.h`) number of characters.

#### Return Value

The function returns `ERR_Ok` in case of success or `ERR_UnknownErrorCode` if the error code is unknown.

#### Remarks

It provides an ASCII string with the description of each error code.

#### Required Header

`PhCon.h`

#### 5.12.3. PhSetDllsLogOption

```
HRESULT PhSetDllsLogOption(UINT SelectCode, UINT Value);
```

#### Parameters

*SelectCode* [in]

It selects the module to set the log option for.

To select a certain module or option use the constants defined in `PhCon.h`:

Value	Meaning
<code>SDLO_PHANTOM</code>	Phantom application logging
<code>SDLO_PHCON</code>	PhCon module logging
<code>SDLO_PHINT</code>	PhInt module logging
<code>SDLO_PHFILE</code>	PhFile module logging
<code>SDLO_PHSIG</code>	PhSig module logging
<code>SDLO_PHSIGV</code>	PhSigV module logging
<code>SDLO_TORAM</code>	Logging to RAM. Faster, but lost if application crashes.

*Value* [in]

Value	Meaning
<code>0xFFFFFFFF</code>	Enable all log dumps from a module
<code>0</code>	Disable all log dumps

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function allows enabling and disabling logging for each module. For example, to enable logging for the `PhCon` module, use the `SDLO_PHCON` selector with *Value* set to `0xFFFFFFFF`. To stop logging use the same selector with *Value* set to `0`.

These settings are also accessible in Help/About/Logging from Phantom or in Manager/Preferences from PCC.

These settings are stored in Windows registry and are unique for all applications using the Phantom SDK.

For a minimum impact on application performance, set `SDLO_TORAM` to `1`. However, please note that if the application crashes or fails to unregister, log information may be lost. To debug crashes set `SDLO_TORAM` to `0`.

#### Required Header

`PhCon.h`

#### 5.12.4. PhGetDllsLogOption

```
HRESULT PhGetDllsLogOption(UINT SelectCode, PUINT pValue);
```

#### Parameters

*SelectCode*

It selects the module to get log options for. See `PhSetDllsLogOption` for possible values.

*pValue*

Pointer where to store requested log option value.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Call this function to find the current logging settings. These settings are used by all applications build on Phantom libraries and are persistent.

#### Required Header

`PhCon.h`

### 5.13. Other - Image Handling & Processing

#### 5.13.1. PhiHtoBITMAPINFOHEADER

```
HRESULT PhiHtoBITMAPINFOHEADER(PIH pIH, PBITMAPINFOHEADER pBMIH);
```

#### Parameters

*pIH* [in]

Buffer containing input image header to be converted.

*pBMIH* [out]

Buffer where to store the `BITMAPINFOHEADER` equivalent of the input image header.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

This function converts a `IH` structure to a `BITMAPINFOHEADER` structure. The two input pointers may be identical in which case the conversion is done in place. Use this function whenever a `BITMAPINFOHEADER` structure is needed for describing one or more images for which only a `IH` structure is available. Please note that even if `IH` is an extension of `BITMAPINFOHEADER` structure it is not safe to use a cast for converting from one type to the other.

### Required Header

`PhInt.h`

### 5.13.2. PhProcessImage

```
HRESULT PhProcessImage(
    PBYTE pPixelSrc,
    PBYTE pPixelDest,
    PIH pIH,
    UINT ProcID,
    PVOID pProcParams
);
```

### Parameters

*pPixelSrc* [in]

Input pixel data to be processed.

*pPixelDest* [out]

Output buffer containing processed data.

*pIH* [in, out]

Input: image header describing input data.

Output: image header describing output data.

*ProcID* [in]

Numerical value indicating the processing operation to be done. Possible values are:

Value	Meaning
<code>IMG_PROC_REDUCE16TO8</code>	Reduce image bit depth to 8 bits per sample

*pProcParams* [in]

Pointer to image processing parameters. The image processing parameters can be of one of the following data types:

<i>procID</i>	Data type
<code>IMG_PROC_REDUCE16TO8</code>	<code>REDUCE16TO8PARAMS</code>

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

The function processes the input cine as indicating by the numerical selector *ProcID* and according to *pProcParams* parameters. *pPixelSrc* and *pPixelDest* may be identical, in which case, the operation is done in place.

Use `PhIntToBitmapInfoHeader` function if the image is used in a context where a `BITMAPINFOHEADER` is requested.

## Required Header

`PhInt.h`

### 5.13.3. `PhInterpolateColor`

```
void PhInterpolateColor(
    PBITMAPINFOHEADER pBMIH,
    PBYTE pPixel,
    UINT CFA,
    UINT Algorithm
);
```

#### Parameters

*pBMIH* [in, out]

Input: image header describing input data.

Output: image header describing output data.

*pPixel* [in, out]

Input: Raw image data.

Output: Color interpolated image data.

*CFA*

Color filter array for the sensor on the camera this image was acquired with. Possible values are:

Value	Meaning
CFA_BAYER	GB RG
CFA_BAYERFLIP	RG GB
CFA_VRI	GBRG RGGB
CFA_VRIV6	BGGR GRBG

To obtain a camera CFA type call `PhGetVersion` with the `GV_CFA` selector.

To obtain the CFA for a cine identified based on a `CINEHANDLE` object, call `PhGetCineInfo` with the `GCI_CFA` selector.

#### Algorithm

Numerical value identifying the algorithm to be used for demosaicing the input image. Different speed/quality combinations are available (see `PhInt.h`):

`FAST_ALGORITHM`

`BEST_ALGORITHM`

`NO_DEMOSAICING`

## No Return Value

## Remarks

Builds a color RGB image from an uninterpolated input image containing raw sensor pixels. Both 8 bits per sample and 16 bits per sample input images can be processed by this function. After demosaicing, the *pBMIH* is updated correspondingly.

Since the color interpolation is done in place, the input *pPixel* buffer must be three times larger than the actual input raw data.

When interpolating an image from a multi-head Phantom v6 camera, this function is able to recognize the gray heads and to interpolate only the

color quarters. This is done using the most significant bits of the CFA selected using the TL\_GRAY, TR\_GRAY, BL\_GRAY, BR\_GRAY masks. The pixels from the gray heads will have their gray level copied three times in the result for the Blue, Green and Red components.

## Required Header

PhInt.h

### 5.13.4. PhTriplicateGray

```
EXIMPROC void PhTriplicateGray(
    PBITMAPINFOHEADER pBMIH,
    PBYTE pPixel,
);
```

#### Parameters

*pBMIH*

DIB header.

*pPixel*

Pixel array.

#### No Return Value

#### Remarks

Build a color RGB image from a monochrome bitmap by copying the same value in the three color component. The input image can be 8 or 16 bit per pixel, the result will be 24 or 48 bpp. The biBitCount and biSizeImage fields of the header are updated.

The pPixel array must be allocated for the resulting color image that is three times larger than the original array. This function is needed in LabVIEW.

#### Required header

PhInt.h

### 5.13.5. PhWBAdjustUnint

```
void PhWBAdjustUnint(
    PBITMAPINFOHEADER pBMIH,
    PBYTE pPixel,
    UINT CFA,
    PWBGAIN pWBg
);
```

#### Parameters

*pBMIH* [in]

Image header describing input data.

*pPixel* [in, out]

Input: Raw image data.

Output: White balance corrected raw image data.

*CFA*

Color filter array for the sensor on the camera this image was acquired with. Possible values are:

Value	Meaning
-------	---------



CFA_BAYER	GB RG
CFA_BAYERFLIP	RG GB
CFA_VRI	GBRG RGGB
CFA_VRIV6	BGGR GRBG

To obtain a camera CFA type call `PhGetVersion` with the `GV_CFA` selector.  
To obtain the CFA for a cine identified based on a `CINEHANDLE` object, call `PhGetCineInfo` with the `GCI_CFA` selector.

*pWBg*

Gain corrections on red and blue channels used to adjust the white balance.

### No Return Value

### Remarks

`PhWBAdjustUnint` applies white balance gain corrections for a raw image.  
For better results, white balance correction must be run before color interpolation.

### Required Header

`PhInt.h`

### 5.13.6. PhComputeWB

```
HRESULT PhComputeWB(
    PBITMAPINFOHEADER pBMIH,
    PBYTE pPixel,
    PPOINT pP,
    int SquareSide,
    UINT CFA,
    PWBGAIN pWB,
    PUINT pSatCnt
);
```

### Parameters

*pBMIH* [in]

Image header.

*pPixel* [in]

Image pixels.

*pP* [in]

The point around which white balance is to be computed.

*SquareSide* [in]

Side length of the square area used to compute white balance.

*CFA* [in]

Color filter array for the sensor on the camera this image was acquired with.

*pWB* [out, optional]

Address where to store corrective gains for red and blue.  
This parameter can be `NULL` if not needed.

*pSatCnt* [out, optional]

Address where to store saturated pixel count.  
This parameter can be `NULL` if not needed.

### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

### Remarks

It computes the compensating `WBGAIN` based on the pixels color in the specified square. For the same square, it can also return saturated pixel count.

### Required Header

`PhCon.h`

### 5.13.7. PhImHisto

```
void PhImHisto(
    PBITMAPINFOHEADER pBMIH,
    PBYTE pPixel,
    UINT HistoSize,
    UINT Color,
    PUINT pHisto,
    PUINT pAvgR,
    PUINT pAvgG,
    PUINT pAvgB
);
```

### Parameters

*pBMIH* [in]

Image header.

*pPixel* [in]

Input image data.

*HistoSize* [in]

Number of elements for the input histogram array.

Can be less than the number of possible levels for the selected color; in this case `[0, MaxPix]` is mapped to `[0, HistoSize - 1]`.

*Color* [in]

Compute the histogram on this selected color:

Value	Meaning
0	RGB converted to gray.
1	Red.
2	Green.
3	Blue.

*pHisto* [out]

Array of integers counting the occurrence of each level or set of levels (if *HistoSize* < maximum value of the pixel).

*pAvgR* [out, optional]

Average Red channel value for the entire image.

This parameter can be `NULL` if not needed.

*pAvgG* [out, optional]

Average Green channel value for the entire image.

This parameter can be `NULL` if not needed.

*pAvgB*

Average Blue channel value for the entire image.

This parameter can be `NULL` if not needed.

### No Return Value

#### Remarks

It computes the histogram for the input image.

#### Required Header

PhInt.h

## 5.14. LabView Specific Functions

### 5.14.1. PhLVRegisterClientEx

```
HRESULT PhLVRegisterClientEx(
    char *pStgPath,
    PROGRESSCALLBACK pfnCallback,
    UINT PhConHeaderVersion
);
```

#### Parameters

*pStgPath* [in]

The path for the stg files of cameras.

It may be your exe path or your StartIn path.

*pfnCallback* [in, optional]

Pointer to a caller function that will be called periodically to indicate the progress during the search for the connected cameras.

This pointer can be NULL if a progress indicator is not needed or it can be an integer from 1 to 15 if PhLVProgress use instead of a callback is preferred.

*PhConHeaderVersion* [in]

PhCon library version. Pass the value of PHCONHEADERVERSION symbol defined in PhCon.h.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

This function registers a LabView application as a client of the Phantom dlls. Call `PhLVRegisterClientEx` before any other function (except for `PhSetDllsOption`) whether working with Phantom camera or with cine files.

#### Required Header

PhCon.h

### 5.14.2. PhLVUnregisterClient

```
HRESULT PhLVUnregisterClient();
```

#### Parameters

No parameters.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

## Remarks

The function unregisters a LabView application previously registered as a client of the Phantom DLLs. Call this function to ensure proper resource deallocation when the Phantom DLLs are not needed anymore.

## Required Header

PhCon.h

### 5.14.3. PhLVProgress

```
HRESULT PhLVProgress(
    UINT CallID,
    PUINT pPercent,
    BOOL Continue
);
```

## Parameters

*CallID* [in]

The call identifier (a number from 1 to 15), which was passed as a *PROGRESSCALLBACK* parameter in a time consuming function.

*pPercent* [out]

The percentage of work done for the time consuming routine call, identified by *CallID*.

*Continue* [in]

Specifies if the function must continue.

Value	Meaning
TRUE	Continue
FALSE	Stop the function call, and force it to return

## Return Value

The function returns *ERR\_Ok* in case of success or an error code, otherwise. To get the text of the error message call *PhGetErrorMessage*.

## Remarks

This function is used as a progress indicator for time-consuming functions. It is designed to work with environments which do not support function pointers (like LabView).

The time-consuming function call and the *PhLVProgress* call must be in different threads.

In LabView the time-consuming function call node must be set to Reentrant.

All API calls that have a *PROGRESSCALLBACK* function pointer parameter accept a *CallID* in it. *CallID* is a *UINT* and it must be in the range [1..15].

For example, the functions can be imagined as having the following headers:

```
// Phcon.dll
HRESULT PhBlackReference(UINT CN, UINT CallID);
HRESULT PhNVMErase(UINT CN, UINT CallID);
HRESULT PhNVMRestore(UINT CN, UINT CineNo, UINT CallID);
HRESULT PhNVMSave(UINT CN, INT Cine, UINT CallID);
HRESULT PhNVMSaveClip(UINT CN, int Cine, PIMRANGE pRng, UINT
Options, UINT CallID);

// Phfile.dll
HRESULT PhWriteCineFile(CINEHANDLE hC, UINT CallID);
```

After calling one of these functions, you can call `PhLVProgress` from another thread to see the amount of work done or to stop the operation.

Example:

```
// Thread 1
PhWriteCineFile(hC, 1); //time consuming function
// Thread 2
UINT Percent;
//...
while(true)
{
    //see the progress of PhWriteCineFile
    PhLVProgress(1, &Percent, TRUE);
    Sleep(1000);
}
```

Important:

To make this function work, the application must register itself by calling `PhLVRegisterClientEx`.

#### Required Header

`PhCon.h`

## 5.15. Nucleus Specific Functions

### 5.15.1. PhUpgradeFirmware

```
HRESULT PhUpgradeFirmware(
    UINT IpAddress,
    UINT CamVer,
    UINT Selector,
    char *pFilename);
```

#### Parameters

*IpAddress* [in]

Camera Ip Address

*CamVer* [in]

Camera hardware version

*Selector* [in]

Selector value indicating the firmware to be upgraded.

*pFilename* [in]

Pointer to complete file name and path. The file must be a .Phfw.

#### Return Value

The function returns `ERR_Ok` in case of success or an error code, otherwise. To get the text of the error message call `PhGetErrorMessage`.

#### Remarks

Use this function to upgrade a camera's firmware, FPGA, (Field Programmable Gate Array), Flash FPGA, Cinemag, and Kernel firmware. For a complete list of currently available selectors, go to Appendix 7.

#### Note:

- 1) The external third-party applications "plink.exe" and "pscp.exe" are required for this function. They should be in the same folder as the application and PhCon.Dll.
- 2) After a camera upgrade, an application program restart should be performed to resynchronize any application cache settings with the updated camera.

#### Required Header

PhCon.h

## 6. File Naming Convention

### The "language" for specifying file names

When naming files saved from the Phantom software, you can use a few special seldom used characters that are legal for filenames. Phantom software recognizes these wildcard characters and substitutes them with different characters in order to better support batch processing.

#### @ followed by a single digit (1 to 8)

Used to insert the cine file number to the path and/or file names. The digit indicates the number of digit counter places the numbering scheme should use. The @ character can be used within both path names and file names. The cine number counter is initialized to zero each time the software is started and incremented by 1 each time a cine is saved to a drive. Use it when running Continuous Recording.

#### + followed by a single digit (1 to 8)

Used to insert an image number in the image file name, starting from 1. The digit number specifies how many digit counter places will be used for the image number.

#### ! followed by a single digit (1 to 8)

Used to insert an image number in the file name, starting from the first frame's number. The digit number specifies how many digit counter places will be used for the image number. Remember to consider the minus sign in the image number, too.

#### \$ followed by a single digit (1 to 8)

Used to denote a session number, if desired, in the path and/or file names. The digit number indicates how many digit counter places the numbering scheme should use.

#### ~ followed by a single digit (8 only)

Used to insert the current date and time within the file names only. The digit number "8" specifies the number of digit places used to report date in the file.

In Continuous Recording character ~ can also be used, like the following:

- a ~ character followed by anything else but 8 will add the cine full trigger time string, including day of week, year, month, .... second, fractions.
- the ~8 characters after the root file name will add 8 characters that represent part of the cine trigger time string: first two characters from the month name, the day, the hour and the minute.

### 6.1. Converting Cines into Multiple Images

You can convert one or several cines into multiple images using this "language", according to your needs. The destination image files can be stored in different ways.

#### 6.1.1. Converting Single Cine into Multiple Images

- When you are converting only one cine, the image files can be all saved in the same folder where the source cine is placed. In this case, you have to

mention just the desired name of the file, followed by a + or a ! and a number. The file names will count from 1 or, respectively, from the first frame's number.

- If you want the image files to be saved in a new folder or a series of new folders, you have to mention the name of the new folders followed by a \ and the name of the image files (described above) and the Phantom software will create them as subfolders within the one in which the source cine is placed.
- The name of the image files described above can also be preceded by a full path name, containing one or several new folders. Phantom software will create the full path before saving the image files.

### 6.1.2. Converting Multiple Cines into Multiple Images

The same rules apply in converting several cines, too, but in this case the Phantom software creates new folders for each cine, named like the source cines (except the file extension), each one containing its own destination image files. These new folders are subfolders in the one in which the source cines are placed or they are placed in new folders according to the complete name given by you.

The names of your new folders can also contain a counter of cines – by using the @ character – or a counter of sessions – by using the \$ character – or the date-and-time string – by using the ~ character.

## 6.2. Examples

1. When converting a cine into multiple images – let's say `cine1.cin` – containing a range of frames from -10 to 2, you can specify the generic file name like this:

```
imag!4
```

The images will be saved in the same folder with the one in which `cine1.cin` is, in files named:

```
imag-010
imag-009
imag-008
...
imag0000
imag0001
imag0002
```

2. Besides the file names, you can indicate a new directory or the full path name, using the wildcards described above. For example, when converting several cines into multiple images – `cine1.cin`, `avi1.avi` – containing a range of frames from -10 to 2, you can specify the following Files in the Multiple Convert Destination dialog box:

```
c:\cines\imag!4
```

The Phantom software will create `cines` folder if it doesn't exist, it will insert new folders named the same like the source files (without extension) – `cine1` and `avi1` in this example – and the images will be saved like this:

```
c:\cines\cine1\imag-010
c:\cines\cine1\imag-009
...
```



```
c:\cines\cinel\imag0001
c:\cines\cinel\imag0002

c:\cines\avi1\imag-010
c:\cines\avi1\imag-009
...
c:\cines\avi1\imag0001
c:\cines\avi1\imag0002
```

3. You can also introduce a counter of cines, using the @ sign. For the same conversion as the one mentioned above, if you specify:

```
c:\cine@2\imag!4
```

you will get the following:

```
c:\cine00\cinel\imag-010
c:\cine00\cinel\imag-009
...
c:\cine00\cinel\imag0001
c:\cine00\cinel\imag0002

c:\cine01\avi1\imag-010
c:\cine01\avi1\imag-009
...
c:\cine01\avi1\imag0001
c:\cine01\avi1\imag0002
```

4. If you prefer a counting of images starting from 1, here is an example for converting the file cinel.cin with the frame range from -10 to 2:

```
c:\cines\imag+4
```

and you will get:

```
c:\cines\imag0001
c:\cines\imag0002
...
c:\cines\imag00012
c:\cines\imag00013
```

5. Continuous recording file naming can be simply defined like the following:

```
cont@4
```

and the cines will be saved into:

```
cont0001
cont0002
cont0003
...
```

If you wish to add the trigger time string to the file name, you can define it:

```
cont ~
```

and you will get:

```
cont Mon Dec 06 2010 16 51 19.274 658
cont Mon Dec 06 2010 16 51 26.009 567
...
```

Phantom software adds WeekDay Month Day Year Hour Minute Second Fractions of second to the file name root.

If you define the name like this:

```
cont ~8
```

the files will be named:

cont DE061710

cont DE061711

...

containing the first two characters of the month name, the day, the hour and the minute.

## 7. Supported File Formats

### Supported cine formats:

FileType	Code	Value	Extension	Write	Read
Cine Raw	MIFILE_RAWCINE	0	.cine	✓	✓
Cine	MIFILE_CINE	1	.cin	✓	✓
AVI	MIFILE_AVI	3	.avi	✓	✓
Multipage TIFF	MIFILE_TIFCINE	4	.tif	✓	✗
QTIME (uncompressed)	MIFILE_QTIME	8	.mov	✓	✗

### Supported image formats:

FileType	Code	Value	Extension	Write	Read
TIFF 8/24 bpp	SIFILE_TIF8	-8	.tif	✓	✓
TIFF 12/36 bpp	SIFILE_TIF12	-9	.tif	✓	✓
TIFF 16/48 bpp	SIFILE_TIF16	-10	.tif	✓	✓
RAW	SIFILE_RAW	-25	.raw	✓	✗
DNG	SIFILE_DNG	-26	.dng	✓	✗
DPX	SIFILE_DPX	-27	.dpx	✓	✗

## 8. Error and Warning Codes

### 8.1. Error Codes

Value	Label	Meaning
-200	ERR_NULLPointer	A NULL pointer was provided for an input parameter or for an output parameter that must be non-NULL.
-201	ERR_MemoryAllocation	Memory allocation failed.
-202	ERR_NoWindow	The hWnd parameter is not a correct window handle.
-203	ERR_CannotRegisterClient	Cannot register client.
-204	ERR_CannotUnregisterClient	Cannot unregister client.
-205	ERR_AsyncRead	1394 communication error.
-206	ERR_AsyncWrite	1394 communication error.
-207	ERR_IsochCIPHeader	1394 isochronous transfer error (Common Isochronous Packet header error).
-208	ERR_IsochDBCContinuity	1394 isochronous transfer error (Data Block Continuity error).
-209	ERR_IsochNoHeader	1394 isochronous transfer error (Sony Vaio).
-210	ERR_IsochAllocateResources	1394 isochronous transfer error: resource allocation.
-211	ERR_IsochAttachBuffers	1394 isochronous transfer error: attach buffers.
-212	ERR_IsochFreeResources	1394 isochronous transfer error: free resources.
-213	ERR_IsochGetResult	1394 isochronous transfer error: get result.
-214	ERR_CannotReadTheSerialNumber	The serial number of the camera cannot be read. The serial number is an important element in the management of the 1394 bus so this is a serious error.
-215	ERR_SerialNumberOutOfRange	The serial read from the camera is incorrect and should be restored.
-216	ERR_UnknownCameraVersion	The connected hardware is not supported by this version of PhCon.
-217	ERR_GetImageTimeout	The image from camera was not available in the allowed time.
-218	ERR_ImageNoOutOfRange	A recorded image has the maximum number PTFrames-1 and the minimum number PTFrames-ImCount.
-220	ERR_CannotReadStgHeader	Stg file read error or incorrect content.
-221	ERR_ReadStg	Stg file read error or incorrect content.
-222	ERR_StgContents	Stg file read error or incorrect content.
-223	ERR_ReadStgOffsets	Stg file read error or incorrect content.
-224	ERR_ReadStgGains	Stg file read error or incorrect content.

Value	Label	Meaning
-225	ERR_NotAStgFile	Stg file read error or incorrect content.
-226	ERR_StgSetupChecksum	Stg file read error or incorrect content.
-227	ERR_StgSetup	Stg file read error or incorrect content.
-228	ERR_StgHardAdjChecksum	Stg file read error or incorrect content.
-229	ERR_StgHardAdj	Stg file read error or incorrect content.
-230	ERR_StgDifferentSerials	The Stg file name is different from the serial stored inside.
-231	ERR_WriteStg	Write to Stg file error.
-232	ERR_NoCine	A recorded image, block or a time buffer was requested but this camera didn't record a cine.
-233	ERR_CannotOpenDevice	1394 isochronous transfer error: cannot open the device.
-234	ERR_TimeBufferSize	The time buffer reported by the camera has an unexpected size.
-236	ERR_CannotWriteCineParams	Cannot write Cine Params to the single cine camera (v4-v6).
-250	ERR_NVMEError	Non volatile memory error.
-251	ERR_NoNVM	Non volatile memory not installed.
-253	ERR_FlashEraseTimeout	Flash erase timeout.
-254	ERR_FlashWriteTimeout	Flash write timeout.
-255	ERR_FlashContents	Flash content is bad.
-256	ERR_FlashOffsetsChecksum	Flash offsets checksum bad.
-257	ERR_FlashGainsChecksum	Flash gains checksum bad.
-258	ERR_TooManyCameras	Too many cameras connected.
-259	ERR_NoResponseFromCamera	No response from Ethernet camera.
-260	ERR_MessageFromCamera	Camera produced an error message.
-261	ERR_BadImgResponse	Bad response to the image request.
-262	ERR_AllPixelsBad	Calibration: all sensor pixels seem to be bad.
-263	ERR_BadTimeResponse	Bad response to the time request.
-264	ERR_GetTimeTimeOut	Time out while getting the image time buffer.
-270	ERR_InBlockTooBig	BASE64 transfer: input block too big.
-271	ERR_OutBufferTooSmall	BASE64 transfer: buffer too small.
-272	ERR_BlockNotValid	BASE64 transfer: Block not valid.
-273	ERR_DataAfterPadding	BASE64 transfer: Data after padding.
-274	ERR_InvalidSlash	BASE64 transfer: Invalid slash.

Value	Label	Meaning
-275	ERR_UnknownChar	BASE64 transfer: Unknown character.
-276	ERR_MalformedLine	BASE64 transfer: Malformed line.
-277	ERR_EndMarkerNotFound	BASE64 transfer: End marker not found.
-280	ERR_NoTimeData	No time recorded.
-281	ERR_NoExposureData	No exposure recorded.
-282	ERR_NoRangeData	No range data recorded.
-283	ERR_NotIncreasingTime	The image time does not increase from one image to the next.
-284	ERR_BadTriggerTime	The trigger time is bad.
-285	ERR_TimeOut	Timeout at the communication to the camera.
-286	ERR_NullWeightsSum	The sum of all weights for partitions size cannot be 0.
-287	ERR_BadCount	Count of partitions cannot be 0 or > maximum partition count.
-288	ERR_CannotChangeRecordedCine	Cannot change the parameters of a recorded cine.
-289	ERR_BadSliceCount	Count of visible i3 slices is different from the stored value.
-290	ERR_NotAvailable	The device you requested information about is not present in the system.
-291	ERR_BadImageInterval	The interval between images is not constant and equal to 1/FrameRate.
-292	ERR_BadCameraNumber	The camera number cannot be above the camera count.
-293	ERR_BadCineNumber	The cine number parameter is incorrect.
-294	ERR_BadSyncObject	A mutex, semaphore or other sync object cannot be created.
-295	ERR_IcmpEchoError	Icmp (ping) echo error.
-296	ERR_MlairReadFirstPacket	Error reading first packet from file.
-297	ERR_MlairReadPacket	Error reading packet from file.
-298	ERR_MlairIncorrectOrder	Incorrect ordered packets detected.
-299	ERR_MlairStartRecorder	Error starting MLAIR recorder process.
-300	ERR_MlairStopRecorder	Error stopping MLAIR recorder process.
-301	ERR_MlairOpenFile	Error opening input data file.
-302	ERR_CmdMutexTimeOut	Command mutex time out.
-303	ERR_CmdMutexAbandoned	Mutex abandoned.
-304	ERR_UnsupportedConversion	Unsupported data conversion.
-305	ERR_TenGigLostPacket	10 Giga lost packet.
-306	ERR_TooManyImgReq	Too many simultaneous image requests.
-307	ERR_BadImRange	Bad image range.

Value	Label	Meaning
-308	ERR_ImgBufferTooSmall	Image buffer too small.
-309	ERR_ImgSize0	Image size is 0.
-310	ERR_IppError	Ipp library error.
-2000	ERR_BadCine	Cine handle has not been initialized.
-2002	ERR_UnsupportedFormat	File format unsupported.
-2003	ERR_InsufficientAlloc	The allocated size of the data is too small.
-2004	ERR_NotInRange	Image number is not in the cine image range.
-2005	ERR_SaveAvi	Error in save AVI process.
-2006	ERR_Encoder	The selected encoder cannot be added.
-2007	ERR_UnsupportedTiffFormat	Tiff file format unsupported.
-2009	ERR_SplitQuartersUnsupported	Split Quarters option not available on selected cine format (AVI and multipage Tiff).
-2010	ERR_NotACineFile	The file is not a cine.
-2012	ERR_FileOpen	Cannot open file.
-2013	ERR_FileRead	Cannot read file.
-2014	ERR_FileWrite	Cannot write file.
-2015	ERR_FileSeek	Cannot set file pointer.
-2016	ERR-DecompressImage	Cannot decompress image.
-2017	ERR_CineVerNewer	Cine file version is newer than the supported cine version.
-2018	ERR_NotSupported	File format not supported.
-2021	ERR_CannotBuildGraph	Cannot build filter graph.
-2023	ERR_AuxDataNotFound	Auxiliary Data not found.
-2024	ERR_NotEnoughMemory	Not enough memory to perform the operation.
-2025	ERR_MpegReaderNotFound	MPEG Reader Interface not found.
-2028	ERR_FunctionNotFound	Function not found in a dynamic link library.
-2029	ERR_CineInUse	Cine is already in use for a similar operation.
-2030	ERR_SaveCineBufferFull	Too many write cine operation are already in progress.
-2031	ERR_NoCineSaveInProgress	No save/write operation is currently in progress for this cine.

## 8.2. Warning Codes

Value	Label	Meaning
0	ERR_Ok	No error, the function completed successfully.
100	ERR_SimulatedCamera	There are no cameras connected or the accessed camera number is above the available cameras. This is intended for testing or using the software without cameras connected.
101	ERR_UnknownErrorCode	The function PhGetErrorMessage cannot provide a message for the specified code.
102	ERR_BadResolution	The resolution specified was incorrect and was adjusted to the maximum camera resolution.
103	ERR_BadFrameRate	The frame rate specified was incorrect and was adjusted to the nearest acceptable value (minimum frame rate if it was less than minimum and maximum frame rate if it was greater than maximum).
104	ERR_BadPostTriggerFrames	Post-trigger frames cannot be 0. It is forced to 1.
105	ERR_BadExposure	The exposure duration specified was incorrect and was adjusted to the nearest acceptable value (minimum exposure if it was less than minimum, and maximum exposure if it was greater than maximum).
106	ERR_BadEDRExposure	EDRExposure cannot be larger than Exposure. It is forced to the value of Exposure.
107	ERR_BufferTooSmall	The buffer provided is too small to fit all the camera resolutions.
108	ERR_CannotSetTime	The camera is synchronized to the IRIG signal and the clock cannot be set manually.
109	ERR_SerialNotFound	PhFindCameraNumber cannot find a camera having the specified serial.
110	ERR_CannotOpenStgFile	The serial tag file (containing settings and calibration) does not exist at the path specified in PhRegisterClient function.
111	ERR_UserInterrupt	User interrupt.
112	ERR_NoSimulatedImageFile	No simulated image file.
113	ERR_SimulatedImageNot24bpp	Simulated image file has to have 24 bpp.
114	ERR_BadParam	Function called with improper argument.
115	ERR_FlashCalibrationNewer	The calibration from camera flash is newer than the stg file.
117	ERR_ConnectedHeadsChanged	The heads connected to this multihead camera changed.
118	ERR_NoHead	Selected head is not installed.
119	ERR_NVMNotInstalled	The Non Volatile Memory is not installed in this camera.
120	ERR_HeadNotAvailable	A cine was stored in a v6.2 camera flash and the recording head is not available at this moment. The images will not be corrected.
121	ERR_FunctionNotAvailable	This function is not available or not implemented in this camera or firmware version.



Value	Label	Meaning
122	ERR_Ph1394dllNotFound	Phantom Firewire driver was not installed - the Firewire cameras will not be visible.
123	ERR_oldtNotFound	Dll's for the signal acquisition (oldaapi32.dll, olm32.dll) are not available. Signal acquisition board will not be visible.
124	ERR_BadFRPSteps	The FRPSteps field value is too large.
125	ERR_BadFRPImgNr	FRPImgNr value is not correct. It has to be between 0 and PTFrames-1.
126	ERR_BadAutoExpLevel	The autoexposure level is not correct. It has to be in the interval [0, 255], regardless of the recording bit depth.
127	ERR_BadAutoExpRect	The autoexposure rectangle is not correct.
128	ERR_BadDecimation	The Decimation factor is not in the interval [1, 16].
129	ERR_BadCineParams	The cine parameters read from camera have incorrect values.
130	ERR_IcmpNotAvailable	Icmp library (or function from this dll) not available.
131	ERR_CorrectResetLine	Special correction for the reset line was requested.
132	ERR_CSRDoneInCamera	Current session Reference was done in camera.
133	ERR_ParamsChanged	Camera parameters were changed from another connection.
134	ERR_ParamReadOnly	Parameter is read-only.
135	ERR_ParamWriteOnly	Parameter is write-only.
136	ERR_ParamNotSupported	Parameter not supported.
137	ERR_IppWarning	Ipp library warning.
140	ERR_BadSensorMode	Invalid sensor mode
141	ERR_CameraOffline	Camera is offline
142	ERR_InvalidFile	Invalid File
2000	ERR_8BitDestinationCine	The destination cine format does not support 16 bit images. The result will be on 8 bit.
2001	ERR_VoidRange	Your selected range does not intersect with the cine image range. The save range contains 0 images.
2002	ERR_ReducedRange	Your selected range intersects partially with the cine image range. The save range is reduced.
2003	ERR_Seeking	This file does not support frame seeking.
2004	ERR_UserAbort	You interrupted this operation.
2005	ERR_FreeLibrary	Cannot free a dynamic loaded library.
2006	ERR_LeadToolsNotFound	LeadTools dlls not found. CCI cine format and some image formats are not available.
2007	ERR_EVSNotFound	EVS library for the MXF file format not found. MXF format will not be available.
2008	ERR_BadResampleParam	Bad resample parameter.

Value	Label	Meaning
2009	ERR_BadCropParam	Bad crop parameter.

## Appendix 1. PhGetCineInfo/PhSetCineInfo Selectors

Name	Meaning	Type	Get	Set
<b>General Cine Parameters</b>				
GCI_ISFILECINE	TRUE if this is a file cine.	BOOL	✓	✗
GCI_FROMFILETYPE	File type (see Supported File Formats).	INT	✓	✗
GCI_COMPRESSION	Compression type. Possible values are: 0 – Gray cine 1 – JPEG compressed cine (*.cci) 2 – Uninterpolated cine (RAW)	UINT	✓	✗
GCI_IMAGECOUNT	Saved image count.	UINT	✓	✗
GCI_TOTALIMAGECOUNT	Total recorded image count.	UINT	✓	✗
GCI_FIRSTIMAGENO	First image number for the selected/saved range.	INT	✓	✗
GCI_FIRSTMOVIEIMAGE	First recorded image number.	INT	✓	✗
GCI_MAXIMGSIZE	Maximum size in bytes required for reading and processing one image of this cine.	UINT	✓	✗
<b>Acquisition Parameters</b>				
GCI_CAMERAVERSION	An integer describing camera version.	UINT	✓	✗
GCI_CAMERASERIAL	The serial of the camera that acquired the cine.	UINT	✓	✗
GCI_FRAMERATE	Cine frame rate.	UINT	✓	✗
GCI_EXPOSURENS	Image exposure in nanoseconds.	UINT	✓	✗
GCI_CFA	CFA of the camera this cine was recorded on.	UINT	✓	✗
GCI_IMWIDTH	Current image width. Can be different from image width at acquisition time if a crop or border operation has been applied for the images of this cine.	UINT	✓	✗
GCI_IMHEIGHT	Current image height. Can be different from image height at acquisition time if a crop or border operation has been applied for the images of this cine.	UINT	✓	✗
GCI_IMWIDTHACQ	Image width at acquisition time.	UINT	✓	✗
GCI_IMHEIGHTACQ	Image height at acquisition time.	UINT	✓	✗

Name	Meaning	Type	Get	Set
GCI_LENSDESCRIPTION	Null terminated string containing various lens information such as producer, model or focal range.	char [ MAXSTDSTRSZ ]	✓	✗
GCI_LENSAPERTURE	Aperture f number.	FLOAT	✓	✗
GCI_LENSFOCALLENGTH	Lens focal length (zoom factor).	FLOAT	✓	✗
GCI_RECORDINGTIMEZONE	The time zone active during the recording of the cine.	INT	✓	✗
GCI_TRIGTC	Trigger frame SMPTE time code and user bits represented as a TC value.	TC	✓	✓
GCI_TRIGTCU	Trigger frame SMPTE time code and user bits represented as a TCU value.	TCU	✓	✓
GCI_PBRATE	Video playback rate (fps) active when the cine was captured.	FLOAT	✓	✓
GCI_TCRATE	Playback rate (fps) used for generating SMPTE time code.	FLOAT	✓	✓
GCI_TRIGTIMESEC	The seconds of the trigger time.	UINT	✓	✗
GCI_TRIGTIMEFR	The fractions of the trigger time.	UINT	✓	✗
GCI_POSTTRIGGER	Post-trigger frames.	UINT	✓	✗
GCI_TRIGFRAME	Information about sync image mode: 0 – internal 1 – external 2 – locktoirig.	UINT	✓	✗
GCI_FRAMEDELAYNS	Frame delay in nanoseconds.	UINT	✓	✗
GCI_EDREXPOSURENS	EDR exposure in nanoseconds.	UINT	✓	✗
GCI_AUTOEXPOSURE	Cine auto-exposure.	UINT	✓	✗
GCI_AUTOEXPLEVEL	Level for autoexposure control.	UINT	✓	✗
GCI_AUTOEXPTOP	Top coordinate of the autoexposure rectangle.	UINT	✓	✗
GCI_AUTOEXPLEFT	Left coordinate of the autoexposure rectangle.	UINT	✓	✗
GCI_AUTOEXPBOTTOM	Bottom coordinate of the autoexposure rectangle.	UINT	✓	✗
GCI_AUTOEXPRIGHT	Right coordinate of the autoexposure rectangle.	UINT	✓	✗
GCI_FRPSTEPS	Number of frame rate profile steps.	UINT	✓	✗
GCI_FRPIMGNOARRAY	Array containing image numbers where to change the acquisition frame rate.	INT[ 16 ]	✓	✗
GCI_FRPRATEARRAY	Array containing new values for the acquisition frame rate.	UINT[ 16 ]	✓	✗

Name	Meaning	Type	Get	Set
GCI_FRPSHAPEARRAY	Array containing frame rate evolution shape for each frame rate profile step. Possible values are: 0 – flat 1 – ramp	UINT[16]	✓	✗
<b>Image Format</b>				
GCI_REALBPP	Pixel color depth for the cine.	UINT	✓	✗
GCI_IS16BPPCINE	TRUE if this is a 16 bpp cine.	BOOL	✓	✗
GCI_ISCOLORCINE	TRUE if this is a color cine.	BOOL	✓	✗
<b>Image Processing Parameters</b>				
GCI_WB	White balance gains for raw cines.	WBGAIN	✓	✓
GCI_WBVIEW	White balance gains applied on color interpolated cines.	WBGAIN	✓	✓
GCI_WBISMETA	If TRUE, white balance gains are stored as metadata.	BOOL	✓	✗
GCI_BRIGHT	Image offset. Range of possible values is [-1.0, 1.0]. 0.0 is neutral. 1.0 means adding the maximum value for the current bit depth.	FLOAT	✓	✓
GCI_CONTRAST	Gain image processing. Range of possible values is [0.1, 10.0]. 1.0 is neutral.	FLOAT	✓	✓
GCI_GAMMA	Per-component gamma Range of possible values is [0.1, 10.0]. 1.0 is neutral.	FLOAT	✓	✓
GCI_GAMMAR	Difference between red channel gamma and overall gamma 0.0 is neutral.	FLOAT	✓	✓
GCI_GAMMAB	Difference between blue channel gamma and overall gamma. 0.0 is neutral.	FLOAT	✓	✓
GCI_SATURATION	Saturation Range of possible values is [0.0, 2.0]. 1.0 is neutral.	FLOAT	✓	✓

Name	Meaning	Type	Get	Set
GCI_HUE	Degrees and fractions of degree to rotate the color hue of every image pixel. Range of possible values is [-180.0, 180.0]. 0.0 is neutral.	FLOAT	✓	✓
GCI_FLIPH	Flip the image horizontally.	BOOL	✓	✓
GCI_FLIPV	Flip the image vertically.	BOOL	✓	✓
GCI_ROTATE	Rotate information: 0 No rotation +90 Counterclockwise -90 Clockwise	INT	✓	✓
GCI_FILTERCODE	Image filter code. See Appendix 2 for possible values.	INT	✓	✓
GCI_IMFILTER	User defined convolution filter.	IMFILTER	✓	✓
GCI_INTALGO	Demosaicing algorithm selector. Phantom SDK offers the following options (see PhInt.h): FAST_ALGORITHM BEST_ALGORITHM NO_DEMOSAICING	UINT	✓	✓
GCI_DEMOSAICINGFUNCPTR	Pointer to current demosaicing responsible function.	DEMOSAICINGFUNCPTR	✓	✓
GCI_TONE	Tone description.	TONEDISC	✓	✓
GCI_RESAMPLEACTIVE	If TRUE, resampling is active.	BOOL	✓	✓
GCI_RESAMPLEWIDTH	Desired resample width.	UINT	✓	✓
GCI_RESAMPLEHEIGHT	Desired resample height.	UINT	✓	✓
GCI_CROPACTIVE	If TRUE, crop is active.	BOOL	✓	✓
GCI_CROPRECTANGLE	Rectangle to crop from the input image.	RECT	✓	✓
GCI_GAIN16_8	Sensitivity value used when converting cine images to 8 bits per sample: $\text{pixelOut} = \text{pixelIn} * \text{fGain16\_8} * (2^8 / 2^{\text{bitdepth}})$	FLOAT	✓	✓

Name	Meaning	Type	Get	Set
GCI_PEDESTALR	Red channel offset to be applied after gamma. Range of possible values is [-1.0, 1.0]. 0.0 is neutral. 1.0 means adjusting by the maximum value for the current bit depth.	FLOAT	✓	✓
GCI_PEDESTALG	Green channel offset to be applied after gamma. Range of possible values is [-1.0, 1.0]. 0.0 is neutral. 1.0 means adjusting by the maximum value for the current bit depth.	FLOAT	✓	✓
GCI_PEDESTALB	Blue channel offset to be applied after gamma. Range of possible values is [-1.0, 1.0]. 0.0 is neutral. 1.0 means adjusting by the maximum value for the current bit depth.	FLOAT	✓	✓
GCI_ENABLEMATRICES	TRUE if the user matrix is currently active.	BOOL	✓	✓
GCI_USERMATRIX	User color matrix description.	CMDESC	✓	✓
GCI_FLARE	Offset to be applied before white balance. Range of possible values is [-1.0, 1.0]. 0.0 is neutral. 1.0 means adjusting by the maximum value for the current bit depth.	FLOAT	✓	✓
GCI_CHROMA	Chrominance adjustment applied after gamma. Range of possible values is [0.0, 2.0]. 1.0 is neutral.	FLOAT	✓	✓
<b>Cine Use Case Parameters</b>				
GCI_VFLIPVIEWACTIVE	For UC_VIEW use cases, cine images are vertically flipped if this parameter is TRUE.	BOOL	✓	✓
GCI_UNCALIBRATEDIMAGE	If TRUE, no calibration correction is applied for this cine images. Camera calibration or CSR results will thus be ignored. Bad pixel correction will also be disabled.	BOOL	✓	✓

Name	Meaning	Type	Get	Set
GCI_NOPROCESSING	If TRUE, no image processing is applied for this cine images except for calibration corrections and bad pixel repairing.	BOOL	✓	✓
GCI_BADPIXELREPAIR	If TRUE, bad pixel correction is enabled.	BOOL	✓	✓
<b>Cine Save Options</b>				
GCI_SAVERANGE	Image range to be saved.	IMRANGE	✓	✓
GCI_SAVEFILENAME	Destination file name.	char[MAX_PATH]	✓	✓
GCI_SAVEFILETYPE	Destination file type. See Supported File Formats for possible values.	INT	✓	✓
GCI_SAVE16BIT	If TRUE, save on 16 bits when more than 8 bits are available.	BOOL	✓	✓
GCI_SAVEPACKED	If >0 and if source cine is packed, destination format is also packed. 0 – Not Packed 1 – Packed 10 2 – Packed 12L	UINT	✓	✓
GCI_SAVEXML	If TRUE, an XML file will be created containing destination cine header.	BOOL	✓	✓
GCI_SAVESTAMPTIME	Image time stamp options. Possible values are: 0 – no time stamp 1 – stamp the absolute time 3 – stamp the time from trigger	INT	✓	✓
GCI_SAVEDECIMATION	Decimation value for image save. Value must be an integer value	FLOAT	✓	✓
GCI_SAVEAVIFRAMERATE	Preferred avi playback frame rate	UINT	✓	✓
GCI_SAVEAVICOMPRESSORLIST	String array containing available avi compressors.	char[MAXCOMPRCNT] [MAXSTDSTRSZ]	✓	✗
GCI_SAVEAVICOMPRESSORINDEX	Preferred avi compressor index.	UINT	✓	✓
GCI_SAVEDPXLSB	If TRUE, preferred DPX byte order is least significant byte first (little endian). Otherwise, big endian byte order is used.	BOOL	✓	✓
GCI_SAVEDPXT010BPS	If TRUE, images are converted to 10 bits per sample.	BOOL	✓	✓



Name	Meaning	Type	Get	Set
GCI_SAVEDPXDATAPACKING	Preferred DPX data packing. Possible values are: 0 – Filled to 32-bit method A 1 – Filled to 32-bit method B 2 – Packed into 32-bit words	UINT	✓	✓
GCI_SAVEDPX10BITLOG	If TRUE, images are converted to 10bit log format.	BOOL	✓	✓
GCI_SAVEDPXEXPORTLOGLUT	If TRUE, look-up tables for linear to log and log to linear are saved in separate files.		✓	
GCI_SAVEDPX10BITLOGREFWHITE	White reference for 10bit log format.	UINT	✓	✓
GCI_SAVEDPX10BITLOGREFBLACK	Black reference for 10bit log format.	UINT	✓	✓
GCI_SAVEDPX10BITLOGGAMMA	Gamma for 10bit log format.	FLOAT	✓	✓
GCI_SAVEDPX10BITLOGFILMGAMMA	Film gamma for 10bit log format.	FLOAT	✓	✓
GCI_SAVEQTPLAYBACKRATE	Preferred playback frame rate for MOV QuickTime files.	UINT	✓	✓
GCI_SAVECCIQUALITY	CCI preferred quality in percents. Range of possible values is [1..100].	UINT	✓	✓
<b>Other Parameters</b>				
GCI_CINENAME	Cine name.	char[MAXSTDSTRSZ]	✓	✓
GCI_WRITEERR	Last error encountered during a save operation from this cine.	INT	✓	✗
GCI_TIMEFORMAT	Preferred time format for this cine. Possible values are (see PhFile.h): TF_LT – Local Time TF_UT – Universal Time TF_SMPTE – SMPTE TimeCode	UINT		
<b>Deprecated</b>				
GCI_HEADSERIAL0	Head 0 serial number of the camera that acquired the cine. (Only for v.6.2 cameras).	UINT	✓	✗
GCI_HEADSERIAL1	Head 1 serial number of the camera that acquired the cine. (Only for v.6.2 cameras).	UINT	✓	✗
GCI_HEADSERIAL2	Head 2 serial number of the camera that acquired the cine. (Only for v.6.2 cameras).	UINT	✓	✗

Name	Meaning	Type	Get	Set
GCI_HEADSERIAL3	Head 3 serial number of the camera that acquired the cine. (Only for v.6.2 cameras).	UINT	✓	✗
GCI_ISMULTIHEADCINE	TRUE if this is a multi-head cine.	BOOL	✓	✗
GCI_EXPOSURE	Image exposure in microseconds.	UINT	✓	✗
GCI_FRP1X	The image number where to change the rate for the first frame rate profile step.	INT	✓	✗
GCI_FRP1Y	New frame rate value used for the first frame rate profile step.	UINT	✓	✗
GCI_FRP2X	The image number where to change the rate for the second frame rate profile step.	INT	✓	✗
GCI_FRP2Y	New frame rate value used for the second frame rate profile step.	UINT	✓	✗
GCI_FRP3X	The image number where to change the rate for the third frame rate profile step.	INT	✓	✗
GCI_FRP3Y	New frame rate value used for the third frame rate profile step.	UINT	✓	✗
GCI_FRP4X	The image number where to change the rate for the forth frame rate profile step.	INT	✓	✗
GCI_FRP4Y	New frame rate value used for the fourth frame rate profile step.	UINT	✓	✗

## Appendix 2. Image Filtering

Predefined filter codes and kernels (see PhInt.h):

Value	Code	Kernel
1	PREWITT_3x3_V	$\begin{bmatrix} -1 & 0 & 1 \\ -1 & 0 & 1 \\ -1 & 0 & 1 \end{bmatrix}$
2	PREWITT_3x3_H	$\begin{bmatrix} 1 & 1 & 1 \\ 0 & 0 & 0 \\ -1 & -1 & -1 \end{bmatrix}$
3	SOBEL_3x3_V	$\begin{bmatrix} -1 & 0 & 1 \\ -2 & 0 & 2 \\ -1 & 0 & 1 \end{bmatrix}$
4	SOBEL_3x3_H	$\begin{bmatrix} 1 & 2 & 1 \\ 0 & 0 & 0 \\ -1 & -2 & -1 \end{bmatrix}$
5	LAPLACIAN_3x3	$\begin{bmatrix} -1 & -1 & -1 \\ -1 & 8 & -1 \\ -1 & -1 & -1 \end{bmatrix}$
6	LAPLACIAN_5x5	$\begin{bmatrix} -1 & -3 & -4 & -3 & -1 \\ -3 & 0 & 6 & 0 & -3 \\ -4 & 6 & 20 & 6 & -4 \\ -3 & 0 & 6 & 0 & -3 \\ -1 & -3 & -4 & -3 & -1 \end{bmatrix}$
7	GAUSSIAN_3x3	$\frac{1}{16} * \begin{bmatrix} 1 & 2 & 1 \\ 2 & 4 & 2 \\ 1 & 2 & 1 \end{bmatrix}$
8	GAUSSIAN_5x5	$\frac{1}{571} * \begin{bmatrix} 2 & 7 & 12 & 7 & 2 \\ 7 & 31 & 52 & 31 & 7 \\ 12 & 52 & 127 & 52 & 12 \\ 7 & 31 & 52 & 31 & 7 \\ 2 & 7 & 12 & 7 & 2 \end{bmatrix}$
9	HIPASS_3x3	$\begin{bmatrix} -1 & -1 & -1 \\ -1 & 8 & -1 \\ -1 & -1 & -1 \end{bmatrix}$

Value	Code	Kernel
10	HIPASS_5x5	$\begin{bmatrix} -1 & -1 & -1 & -1 & -1 \\ -1 & -1 & -1 & -1 & -1 \\ -1 & -1 & 24 & -1 & -1 \\ -1 & -1 & -1 & -1 & -1 \\ -1 & -1 & -1 & -1 & -1 \end{bmatrix}$
11	SHARPEN_3x3	$\frac{1}{8} * \begin{bmatrix} -1 & -1 & -1 \\ -1 & 16 & -1 \\ -1 & -1 & -1 \end{bmatrix}$
-1	User filter	

## Appendix 3. PhGet/PhSet Selectors

Name	Meaning	Type	Get	Set
<b>Camera Info</b>				
gsModel	Camera model	char[MAXSTDSTRSZ]	✓	✗
gsIPAddress	Camera current Ethernet address	char[MAXIPSTRSZ]	✓	✗
<b>Video out Parameters</b>				
gsVideoPlayCine	Cine number currently played on video out	INT	✓	✗
gsVideoPlaySpeed	Video out play speed in mHz.	INT	✓	✗
gsVideoOutputConfig	Video output type. Possible values are: 0 – Single-feed mode 1 – Single-feed mode with dual-link 444 2 – Dual-feed mode 3 – Dual-feed mode with 444	INT	✓	✓
<b>Temperatures</b>				
gsSensorTemperature	Sensor temperature	UINT	✓	✗
gsCameraTemperature	Camera temperature	INT	✓	✗
gsThermoElectricPower	Amount of power used by the Peltier device	INT	✓	✗
gsSensorTemperatureThreshold	Sensor temperature threshold	INT	✓	✓
gsCameraTemperatureThreshold	Camera temperature threshold	INT	✓	✓
<b>Lens Parameters</b>				
gsLensFocus	Move the focus ring by the given number of incremental units relative to current focus position.	INT	✗	✓
gsLensAperture	fstop value	double	✓	✓
gsLensApertureRange	Focal length range	double[2]	✓	✗
gsLensDescription	Lens description	char[MAXSTDSTRSZ]	✓	✗
gsLensFocusInProgress	Focus in progress	BOOL	✓	✗
gsLensFocusAtLimit	Focus adjustment has hit a mechanical limit	BOOL	✓	✗
<b>Camera Capabilities</b>				

Name	Meaning	Type	Get	Set
gsSupportsInternalBlackRef	Tells if the camera knows to perform a black calibration	BOOL	✓	✗
gsSupportsImageTrig	Image based auto trigger capability	BOOL	✓	✗
gsSupportsCardFlash	Tells if the camera supports a compact flash	BOOL	✓	✗
gsSupportsMagazine	Tells if the camera supports a CineMag	BOOL	✓	✗
gsSupportsHQMode	Tells if HQ mode is supported	BOOL	✓	✗
gsSupportsGenlock	Tells if video genlock is supported	BOOL	✓	✗
gsSupportsEDR	Tells if EDR is supported	BOOL	✓	✗
gsSupportsAutoExposure	Tells if auto-exposure is supported	BOOL	✓	✗
gsSupportsTurbo	Tells if turbo mode is supported	BOOL	✓	✗
gsSupportsBurstMode	Tells if burst mode is supported	BOOL	✓	✗
gsSupportsShutterOff	Tells if shutter off mode is supported	BOOL	✓	✗
gsSupportsDualSDIOutput	Tells if dual SDI output is present	BOOL	✓	✗
gsSupportsRecordingCines	Tells if this device can record cines. (A CineStation can never record cines.)	BOOL	✓	✗
gsSupportsV444	Tells if HD-SDI 4:4:4 output is available	BOOL	✓	✗
gsSupportsInterlacedSensor	Tells if sensor has interlaced channels.	BOOL	✓	✗
gsSupportsRampFRP	Tells if camera supports frame rate profile using the ramp camera field.	BOOL	✓	✗
gsSupportsOffGainCorrections	Tells if camera implements offset gain corrections internally.	BOOL	✓	✗
gsSupportsQuietMode	Tells if camera supports quiet mode (i.e. turning off the fan)	BOOL	✓	✗
gsMaxPartitionCount	Maximum partition count.	UINT	✓	✗
gsSupportsProgIO	Tells if camera supports programmable I/O	BOOL	✓	
gsSupportsSensorModes	Tells if camera supports sensor modes	BOOL	✓	
gsSupportsBattery	Tells if the camera supports a internal battery	BOOL	✓	✗
gsSupportsExpIndex	Tells if the camera supports exposure index	BOOL	✓	
<b>Camera Current Status Information</b>				
gsHasMechanicalShutter	Mechanical shutter present or not	BOOL	✓	✗
gsHasBlackLevel4	Tells if the camera uses a non-zero value for black level.	BOOL	✓	✗

Name	Meaning	Type	Get	Set
gsHasCardFlash	Tells if the camera has a card flash inserted and mounted	BOOL	✓	✗
gsHas10G	Tells if the camera has a 10G connection up and running	BOOL	✓	✗
<b>Camera Current Parameters</b>				
gsMechanicalShutter	Shutter control: on or off	INT	✓	✓
gsImageTrigThreshold	Amount a pixel value must change in order to be counted as an active pixel for auto trigger purposes	INT	✓	✓
gsImageTrigAreaPercentage	Percentage of the area of the auto trigger region that must be active in order for an auto trigger event to be generated	INT	✓	✓
gsImageTrigSpeed	Number of frames between updates of the auto trigger reference memory	INT	✓	✓
gsImageTrigMode	Auto trigger mode: 0 – auto trigger disabled 1 – the camera will both drive the autotrigger and trigger itself when an autotrigger condition is detected. If the autotrigger signal is pulled low by an external device, the camera will be triggered. 2 – the image changes are analyzed, and when an auto trigger condition is detected, the autotrigger signal is pulled low, as in mode 1, but the camera will not trigger itself. An external device pulling low the autotrigger signal will not trigger the camera either.	INT	✓	✓
gsImageTrigRect	Autotrigger region. Top left corner has (0, 0) coordinates.	RECT	✓	✓
gsAutoProgress	Number of frames still to do by the camera auto process.	INT	✓	✗
gsAutoBlackRef	Auto black reference enabled or not (performing a CSR before each new recording)	INT	✓	✓
gsCardFlashError	Last compact flash error	INT	✓	✗
gsGenlock	Video genlock	INT	✓	✓
gsGenlockStatus	Genlock status	char[MAXSTDSTRSZ]	✓	✗
gsTurboMode	Turbo mode support	INT	✓	✓
gsQuiet	Camera fan: on or off	INT	✓	✓
gsClockPeriod	Obtain the camera clock period as a double value (seconds)	double	✓	

Name	Meaning	Type	Get	Set												
gsPortCount	How many ports (fixed & programmable) are available on this camera (0 =indicates none)	UINT	✓													
gsTriggerEdgeAndVoltage	bit 0 rising edge, bit 1 high voltage (6v threshold instead of 1.5v) <table><tr><td>Bit 0 Value</td><td>Meaning</td></tr><tr><td>0</td><td>Trigger on falling edge</td></tr><tr><td>1</td><td>Trigger on rising edge</td></tr></table> <table><tr><td>Bit 1 Value</td><td>Meaning</td></tr><tr><td>0</td><td>Use 1.5 Voltage threshold</td></tr><tr><td>1</td><td>Use high Voltage threshold (6v threshold instead of 1.5v)</td></tr></table>	Bit 0 Value	Meaning	0	Trigger on falling edge	1	Trigger on rising edge	Bit 1 Value	Meaning	0	Use 1.5 Voltage threshold	1	Use high Voltage threshold (6v threshold instead of 1.5v)	UINT	✓	✓
Bit 0 Value	Meaning															
0	Trigger on falling edge															
1	Trigger on rising edge															
Bit 1 Value	Meaning															
0	Use 1.5 Voltage threshold															
1	Use high Voltage threshold (6v threshold instead of 1.5v)															
gsTriggerFilter	Filter constant (seconds)	double	✓	✓												
gsTriggerDelay	Trigger delay in seconds available at camera. Only valid for cameras with programmable I/O capability	double	✓	✓												
gsSensorModesList	Get list of camera Sensor modes	char[MAXIPSTRSZ]	✓													
gsSensorMode	Get/Set Sensor Mode	UINT	✓	✓												
gsBatteryEnable	Enable/Disable battery operation	BOOL	✓	✓												
gsBatteryCaptureEnable	Enable/Disable battery operation during capture	BOOL	✓	✓												
gsBatteryPreviewEnable	Enable/Disable battery operation during preview	BOOL	✓	✓												
gsBatteryWtrRuntime	Time in seconds the camera will run on battery in WTR if a cine is not triggered.	UINT	✓	✓												
gsBatteryVoltage	Battery voltage level	double	✓	✗												



Name	Meaning	Type	Get	Set																
gsBatteryState	Battery control status	UINT	✓	✗																
	<table><tr><th>Value</th><th>Meaning</th></tr><tr><td>0</td><td>battery not present</td></tr><tr><td>1</td><td>battery charging</td></tr><tr><td>2</td><td>battery charging, voltage above 1/2 capacity threshold;</td></tr><tr><td>3</td><td>battery charged;</td></tr><tr><td>4</td><td>battery discharging;</td></tr><tr><td>5</td><td>battery low;</td></tr><tr><td>16</td><td>battery charging fault.</td></tr></table>				Value	Meaning	0	battery not present	1	battery charging	2	battery charging, voltage above 1/2 capacity threshold;	3	battery charged;	4	battery discharging;	5	battery low;	16	battery charging fault.
	Value				Meaning															
	0				battery not present															
	1				battery charging															
	2				battery charging, voltage above 1/2 capacity threshold;															
	3				battery charged;															
	4				battery discharging;															
	5				battery low;															
16	battery charging fault.																			
When the battery is armed, bit 3 is set (8 is added to the values above).																				
gsBatteryArmDelay	Delay, in seconds, from the moment the camera is placed into WTR and the time the battery is armed.	UINT	✓	✓																
gsBatteryPrevRuntime	Time in seconds the camera will run on battery in preview if a cine is not triggered.	UINT	✓	✓																
gsBatteryPwrOffDelay	Time in seconds the battery still supplies power to the camera after it has been disarmed.	UINT	✓	✓																
gsBatteryReadyGate	Set/Get Battery readygate control	BOOL	✓	✓																
gsExpIndex	Set/Get exposure index ISO value	UINT	✓	✓																
gsExpIndexPresets	Get list of camera ISO values	UINT[10]	✓																	
Network Interface Alias																				
gsEthernetDefaultAddress	Primary/Default IP address	char[MAXIPSTRSZ]	✓	✓																
gsEthernetAddress	Secondary IP address DHCP Control: Enable: Set this to the word "auto" Disable: Remove the word "auto"	char[MAXIPSTRSZ]	✓	✓																
	Note: DHCP is only functional on selective cameras																			
gsEthernetMask	Network mask	char[MAXIPSTRSZ]	✓	✓																

Name	Meaning	Type	Get	Set
gsEthernetBroadcast	Broadcast mask	char[MAXIPSTRSZ]	✓	✓
gsEthernetGateway	IP address of default gateway	char[MAXIPSTRSZ]	✓	✓
gsEthernet10GAddress	10G IP address	char[MAXIPSTRSZ]	✓	✓
gsEthernet10GMask	10G Network mask	char[MAXIPSTRSZ]	✓	✓
gsEthernet10GBroadcast	10G Broadcast mask	char[MAXIPSTRSZ]	✓	✓

## Appendix 4. PhCineGet/PhCineSet Selectors

Name	Meaning	Type	Get	Set
<b>Cine Status Parameters</b>				
cgsVideoMarkIn	Start of the videoplay range	INT	✓	✓
cgsVideoMarkOut	End of the videoplay range	INT	✓	✓
cgsIsRecorded	Cine status: recorded or not	BOOL	✓	✗
cgsHqMode	HQ mode status	BOOL	✓	✓
<b>Burst Mode Parameters</b>				
cgsBurstCount	Number of frames in a burst	UINT	✓	✓
cgsBurstPeriod	Interval between two frames in a burst, in nanoseconds	UINT	✓	✓
<b>Lens Information</b>				
cgsLensDescription	Lens description	char[MAXSTDSTRSZ]	✓	✗
cgsLensAperture	Lens aperture f number	float	✓	✗
cgsLensFocalLength	Lens focal length	float	✓	✗
<b>Acquisition Parameters</b>				
cgsShutterOff	Enable/disable maximum exposure for piv mode	BOOL	✓	✓
<b>Time Information</b>				
cgsTriggerTime	Time of trigger	TIME64	✓	✗
cgsTrigTC	Time code of the trigger frame	TC	✓	✗
cgsPbRate	Playback rate (fps) of the video mode active when the cine was captured	float	✓	✗
cgsTcRate	Time code rate (fps) corresponding to the video mode active when the cine was captured	float	✓	✗
<b>Image Processing Parameters</b>				
cgsWb	White balance gains	WBGAIN	✓	✓
cgsBright	Brightness	float	✓	✓
cgsContrast	Contrast	float	✓	✓
cgsGamma	Per component gamma		✓	✓

Name	Meaning	Type	Get	Set
cgsGammaR	Difference between red channel gamma and overall gamma	float	✓	✓
cgsGammaB	Difference between blue channel gamma and overall gamma	float	✓	✓
cgsSaturation	Saturation	float	✓	✓
cgsHue	Hue	float	✓	✓
cgsPedestalR	Red channel pedestal	float	✓	✓
cgsPedestalG	Green channel pedestal	float	✓	✓
cgsPedestalB	Blue channel pedestal	float	✓	✓
cgsFlare	Flare adjustment	float	✓	✓
cgsChroma	Chroma gain	float	✓	✓
cgsTone	A label plus a set of control points describing the tone	TONEDESC	✓	✓
cgsUserMatrix	A label plus a 3 by 3 matrix describing user matrix.	CMDESC	✓	✓
cgsEnableMatrices	Disable/enable user matrix	BOOL	✓	✓
<b>Other Parameters</b>				
cgsName	Cine name	char[MAXSTDSTRSZ]	✓	✓
<b>Deprecated</b>				
cgsVideoTone	Integer selector among predefined video tones	INT	✓	✓

## Appendix 5. PhGet1/PhSet1 Selectors

Name	Meaning	Type	Get	Set
gsSigCount	Signal count available at the port	UINT	✓	✗
gsSigSelect	Get/Select the signal to connect to a certain port	UINT	✓	✓✗
gsPulseProc	Parameters for pulse processor (in a structure)	PULSEPARAM	✓	✓✗
gsHasPulseProc	Tells if this port supports pulse processor control	BOOL	✓	✗

## Appendix 6. PhGet2 Selectors

Name	Meaning	Type	Get	Set
gsSigName	Signal name, given the port number and signal number	Char[ <b>MAXSTDSTRSZ</b> ]	✓	✗

## Appendix 7. PhUpgradeFirmware Selectors

Name	Meaning	Type	Get	Set
FW_PHFW	Upgrade Ph16 camera firmware (.phfw file)	Char[ <b>MAXSTDSTRSZ</b> ]		✓✗